

**Model 672
Spectroscopy Amplifier
Operating and Service Manual**

U.S.A. Patent No. 4,866,400

Advanced Measurement Technology, Inc.

a/k/a/ ORTEC[®], a subsidiary of AMETEK[®], Inc.

WARRANTY

ORTEC* warrants that the items will be delivered free from defects in material or workmanship. ORTEC makes no other warranties, express or implied, and specifically NO WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.

ORTEC's exclusive liability is limited to repairing or replacing at ORTEC's option, items found by ORTEC to be defective in workmanship or materials within one year from the date of delivery. ORTEC's liability on any claim of any kind, including negligence, loss, or damages arising out of, connected with, or from the performance or breach thereof, or from the manufacture, sale, delivery, resale, repair, or use of any item or services covered by this agreement or purchase order, shall in no case exceed the price allocable to the item or service furnished or any part thereof that gives rise to the claim. In the event ORTEC fails to manufacture or deliver items called for in this agreement or purchase order, ORTEC's exclusive liability and buyer's exclusive remedy shall be release of the buyer from the obligation to pay the purchase price. In no event shall ORTEC be liable for special or consequential damages.

Quality Control

Before being approved for shipment, each ORTEC instrument must pass a stringent set of quality control tests designed to expose any flaws in materials or workmanship. Permanent records of these tests are maintained for use in warranty repair and as a source of statistical information for design improvements.

Repair Service

If it becomes necessary to return this instrument for repair, it is essential that Customer Services be contacted in advance of its return so that a Return Authorization Number can be assigned to the unit. Also, ORTEC must be informed, either in writing, by telephone [(865) 482-4411] or by facsimile transmission [(865) 483-2133], of the nature of the fault of the instrument being returned and of the model, serial, and revision ("Rev" on rear panel) numbers. Failure to do so may cause unnecessary delays in getting the unit repaired. The ORTEC standard procedure requires that instruments returned for repair pass the same quality control tests that are used for new-production instruments. Instruments that are returned should be packed so that they will withstand normal transit handling and must be shipped PREPAID via Air Parcel Post or United Parcel Service to the designated ORTEC repair center. The address label and the package should include the Return Authorization Number assigned. Instruments being returned that are damaged in transit due to inadequate packing will be repaired at the sender's expense, and it will be the sender's responsibility to make claim with the shipper. Instruments not in warranty should follow the same procedure and ORTEC will provide a quotation.

Damage in Transit

Shipments should be examined immediately upon receipt for evidence of external or concealed damage. The carrier making delivery should be notified immediately of any such damage, since the carrier is normally liable for damage in shipment. Packing materials, waybills, and other such documentation should be preserved in order to establish claims. After such notification to the carrier, please notify ORTEC of the circumstances so that assistance can be provided in making damage claims and in providing replacement equipment, if necessary.

Copyright © 2002, Advanced Measurement Technology, Inc. All rights reserved.

*ORTEC[®] is a registered trademark of Advanced Measurement Technology, Inc. All other trademarks used herein are the property of their respective owners.

CONTENTS

WARRANTY	ii
SAFETY INSTRUCTIONS AND SYMBOLS	iv
SAFETY WARNINGS AND CLEANING INSTRUCTIONS	v
1. DESCRIPTION	1
1.1. GENERAL	1
2. SPECIFICATIONS	3
2.1. PERFORMANCE	3
2.2. CONTROLS AND INDICATORS	4
2.3. INPUTS	5
2.4. OUTPUTS	5
2.5. ELECTRICAL AND MECHANICAL	6
3. INSTALLATION	7
3.1. POWER CONNECTION	7
3.2. PREAMPLIFIER CONNECTION	7
3.3. PULSED RESET PREAMPLIFIERS AND INHIBIT IN CONNECTION	7
3.4. CONNECTION OF TEST PULSE GENERATOR	7
3.5. SHAPING CONSIDERATIONS	8
3.6. LINEAR OUTPUT CONNECTIONS	8
3.7. PILE-UP REJECTION USING PUR OUTPUT	9
3.8. LIVETIME CORRECTION USING BUSY OUTPUT	9
3.9. INPUT COUNT RATE USING CRM OUTPUT	9
4. OPERATING INSTRUCTIONS	10
4.1. INITIAL TESTING AND OBSERVATION OF PULSE WAVEFORMS	10
4.2. STANDARD SETUP PROCEDURES	10
4.3. POLE-ZERO ADJUSTMENTS FOR RESISTIVE-FEEDBACK PREAMPLIFIER	11
4.4. BASELINE RESTORER (BLR) SETTING	12
4.5. INTERNAL CONTROLS	13
4.6. DIFFERENTIAL INPUT MODE	14
4.7. SYSTEM THROUGHPUT	16
4.8. CHARGE COLLECTION OR BALLISTIC DEFICIT EFFECTS	17
4.9. PILE-UP REJECTOR (PUR) AND LIVETIME CORRECTOR	18
4.10. OPERATION WITH SEMICONDUCTOR DETECTORS	20
4.11. OPERATION IN SPECTROSCOPY SYSTEMS	22
4.12. OTHER EXPERIMENTS	23
5. MAINTENANCE	27
5.1. TEST EQUIPMENT REQUIRED	27
5.2. PULSER TEST	27
5.3. SUGGESTIONS FOR TROUBLESHOOTING	28
5.4. FACTORY REPAIR	28
5.5. TABULATED TEST POINT VOLTAGES	29

SAFETY INSTRUCTIONS AND SYMBOLS

This manual contains up to three levels of safety instructions that must be observed in order to avoid personal injury and/or damage to equipment or other property. These are:

DANGER Indicates a hazard that could result in death or serious bodily harm if the safety instruction is not observed.

WARNING Indicates a hazard that could result in bodily harm if the safety instruction is not observed.

CAUTION Indicates a hazard that could result in property damage if the safety instruction is not observed.

Please read all safety instructions carefully and make sure you understand them fully before attempting to use this product.

In addition, the following symbol may appear on the product:



ATTENTION—Refer to Manual



DANGER—High Voltage

Please read all safety instructions carefully and make sure you understand them fully before attempting to use this product.

SAFETY WARNINGS AND CLEANING INSTRUCTIONS

DANGER Opening the cover of this instrument is likely to expose dangerous voltages. Disconnect the instrument from all voltage sources while it is being opened.

WARNING Using this instrument in a manner not specified by the manufacturer may impair the protection provided by the instrument.


Cleaning Instructions

To clean the instrument exterior:

- Unplug the instrument from the ac power supply.
- Remove loose dust on the outside of the instrument with a lint-free cloth.
- Remove remaining dirt with a lint-free cloth dampened in a general-purpose detergent and water solution. Do not use abrasive cleaners.

CAUTION To prevent moisture inside of the instrument during external cleaning, use only enough liquid to dampen the cloth or applicator.

- Allow the instrument to dry completely before reconnecting it to the power source.

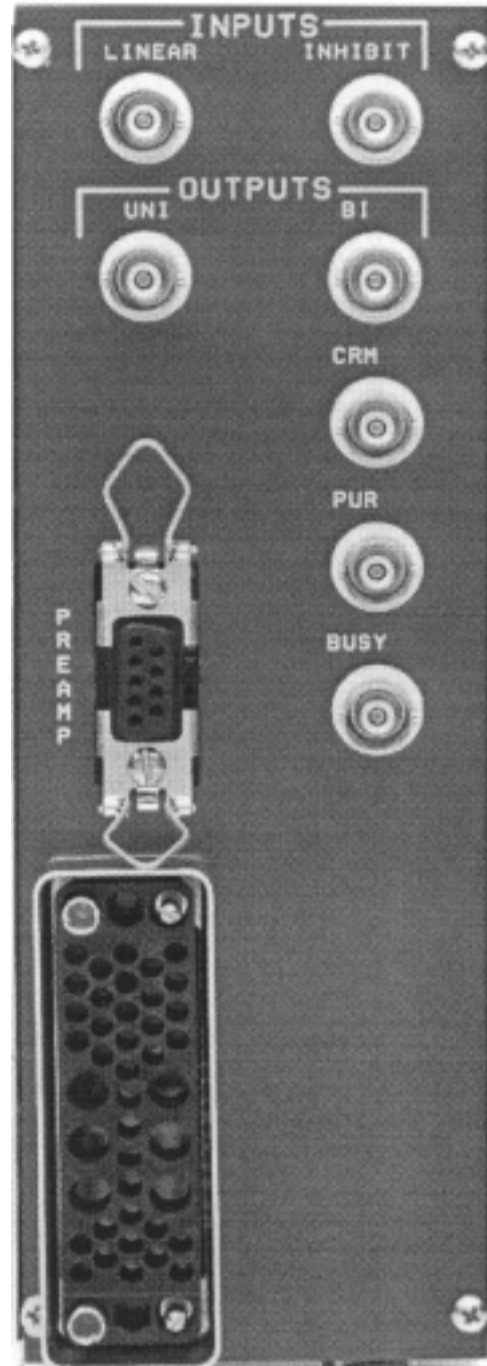

ORTEC MODEL 672 SPECTROSCOPY AMPLIFIER

1. DESCRIPTION

1.1. GENERAL

The ORTEC Model 672 high-performance, energy spectroscopy amplifier is ideally suited for use with germanium, silicon surface-barrier, and Si(Li) detectors. It can also be used with scintillation detectors and proportional counters. The Model 672 input accepts either positive or negative polarity signals from a detector preamplifier and provides a positive 0 to 10-V output signal suitable for use with single- or multichannel pulse height analyzers. Its gain is continuously variable from 2.5 to 1500.

Automation of all the critical adjustments makes the 672 easy to set up with any detector, and provides a performance that is nearly independent of operator expertise.

The Automatic Pole-Zero Adjustment feature significantly simplifies the tuning of the amplifier to compensate for the decay time of the preamplifier pulse. This minimizes the operator skill and effort needed to achieve good energy resolution and peak position stability at moderate to high counting rates (Fig. 1.1.). When changing time constants or detectors, an accurate pole-zero (PZ) adjustment is achieved by simply pushing the AUTO PZ button and waiting a few seconds for the AUTO PZ BUSY LED to turn off. No oscilloscope is required for this procedure. The AUTO PZ memory is protected against unforeseen power failures.

In extreme situations, where the preamplifier pulse shape is deformed from the normal exponential decay, complete PZ cancellation is not possible. In such cases, a slight improvement in the high counting rate performance can sometimes be achieved using the manual PZ adjustment mode to arrive at a compromise solution.

A front-panel switch on the Model 672 provides the choice of either a triangular or a Gaussian pulse shape on the UNIPOLAR output connector. The noise performance of the triangular pulse shape is equivalent to a Gaussian pulse shape having a 17% longer shaping time constant. In applications where the series noise component is dominant (short shaping time constants), and the pile-up rejector is

utilized, the triangular shape will generally offer the same deadtime and slightly lower noise than the Gaussian pulse shape. A front-panel switch permits selection of the optimum shaping time constant for each detector and application. Six time constants in the range of 0.5 to 10 μ s, and the TRIANGULAR/GAUSSIAN switch combine to offer 12 different shaping times. A bipolar output is also provided for measurements requiring zero cross-over timing.

To minimize spectrum distortion at medium and high counting rates (Fig. 2.1.), the unipolar output incorporates a high-performance, gated, baseline restorer with several levels of automation. Automatic positive and negative noise discriminators ensure that the baseline restorer operates only on the true baseline between pulses in spite of changes in the noise level. No operator adjustment of the baseline restorer is needed when changes are made in the gain, the shaping time constant, or the detector characteristics. Negative overload recovery from the reset pulses generated by transistor reset preamplifiers and pulsed optical feedback preamplifiers is also handled automatically. A monitor circuit gates off the baseline restorer and provides a reject signal for a multichannel analyzer until the baseline has safely recovered from the overload.

Several operating modes are selectable for the base-line restorer. For making either a manual or automatic PZ adjustment, the PZ position is selected. This position can also be used where the slowest baseline restorer rate is desired. For situations where low frequency noise interference is a problem, the HIGH rate can be chosen. On detectors where perfect PZ cancellation is impossible, the AUTO baseline restorer rate provides the optimum performance at both low and high counting rates.

An efficient pile-up rejector is incorporated in the 672 Spectroscopy Amplifier. It provides an output logic pulse for the associated multichannel analyzer to suppress the spectral distortion caused by pulses piling up on each other at high counting rates

(Fig. 2.2.). The fast amplifier in the pile-up rejector includes a gated baseline restorer with its own automatic noise discriminator. A multicolor pile-up rejector LED on the front panel indicates the throughput efficiency of the amplifier. At low counting rates the LED flashes green. The LED turns yellow at moderate counting rates and red when pulse pile-up losses are >70%.

When long connecting cables are used between the detector preamplifier output and the amplifier input, noise induced in the cable by the environment can be a problem. The Model 672 provides two solutions. For low to moderate interference

frequencies the differential input mode can be used with paired cables from the preamplifier to suppress the induced noise. At high frequencies a common mode rejection transformer built into the 672 input reduces noise pick-up. The transformer is particularly effective in eliminating interference from the display raster generators in personal computers.

All toggle switches on the front panel lock to prevent accidental changes in the desired settings.

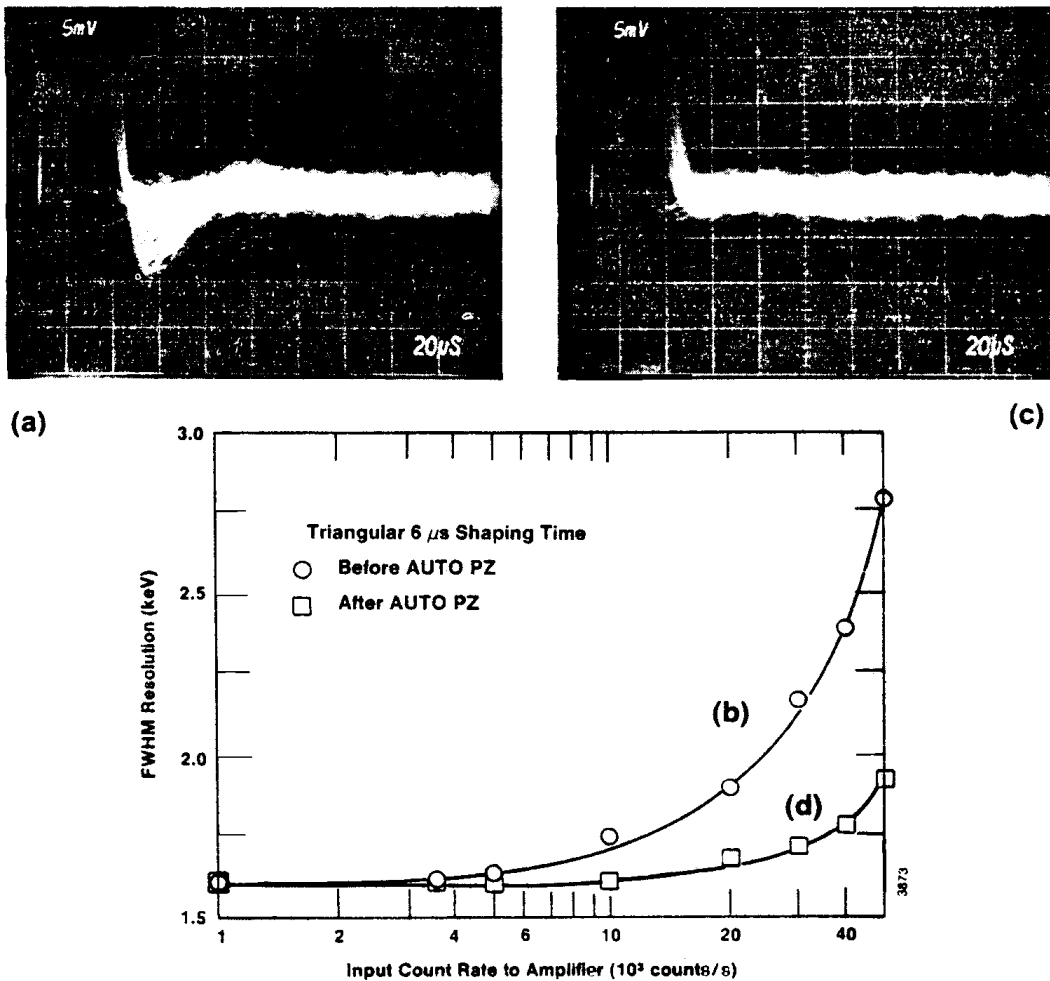


Fig. 1.1. Effectiveness of the Automatic Pole-Zero Feature. (a) An uncanceled pole produces an undershoot on the amplifier output pulse, which, in turn, causes (b) premature degradation of the resolution in the energy spectrum as the counting rate increases. (c) After pushing the AUTO PZ button, the undershoot is automatically removed, resulting in (d) better resolution at high counting rates.

2. SPECIFICATIONS¹

2.1. PERFORMANCE

Note: Unless otherwise stated, performance specifications are measured on the unipolar output with 2- μ s Gaussian shaping, the manual PZ mode, and the AUTO BLR mode.

GAIN RANGE Continuously adjustable from 2.5 to 1500. Gain is the product of the COARSE and FINE GAIN controls.

UNIPOLAR PULSE SHAPES Switch selection of a nearly triangular pulse shape or a nearly Gaussian pulse shape at the UNIPOLAR output (Table 2.1).

BIPOLAR OUTPUT PULSE SHAPE Rise of the BIPOLAR output pulse from 0.1% to maximum amplitude is 1.65 times selected SHAPING TIME. Zero cross-over of the bipolar output pulse delayed from the maximum amplitude of Gaussian UNIPOLAR output by 0.33 times selected SHAPING TIME.

INTEGRAL NONLINEARITY (UNIPOLAR Output)
< $\pm 0.025\%$ from 0 to +10 V.

NOISE Equivalent input noise <5.0 μ V rms for gains >100, and <4.5 μ V rms for gains >1000 in manual PZ mode, or <6.0 μ V for gains >100 in AUTO PZ mode.

TEMPERATURE COEFFICIENT (0 to 50°C)

Unipolar Output < $\pm 0.005\%/^{\circ}\text{C}$ for gain, and < ± 7.5 $\mu\text{V}/^{\circ}\text{C}$ for dc level.

Bipolar Output < $\pm 0.007\%/^{\circ}\text{C}$ for gain, and < ± 30 $\mu\text{V}/^{\circ}\text{C}$ for dc level.

WALK Bipolar zero cross-over walk is < ± 3 ns over a 50:1 dynamic range.

OVERLOAD RECOVERY Unipolar and bipolar outputs recover to within 2% of the rated output from a X 1000 overload in 2.5 non-overloaded pulse widths using maximum gain.

Table 2.1. Unipolar Pulse Shape Parameters for the Triangular and Gaussian Pulse Shapes

Time Interval	Shaping Time Multiplier ^a	
	Triangular	Gaussian
From start of input pulse to maximum amplitude of unipolar output pulse	2.6	2.8
Rise of output pulse from 0.1% to maximum amplitude	2.4	2.0
Width of output pulse at 50% of maximum amplitude	2.5	2.0
Width of output pulse at 1% of maximum amplitude	5.6	5.0
Width of output pulse at 0.1% of maximum amplitude	6.9	6.3

^aTime interval equals the selected front-panel SHAPING TIME multiplied by the Shaping Time Multiplier.

¹ Specifications subject to change without notice.

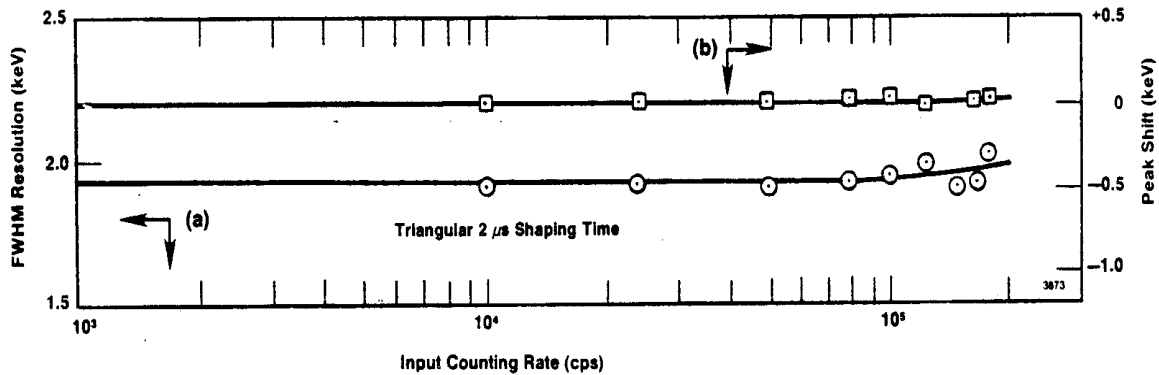


Fig. 2.1. (a) Resolution and (b) Peak Position Stability as a Function of Counting Rate.

SPECTRUM BROADENING† (Fig. 2.1.) Typically <8% broadening of the FWHM for counting rates up to 100,000 counts per second (counts/s), and <15% broadening for counting rates up to 200,000 counts/s. Measured on the 1.33-MeV gamma-ray line from a ⁶⁰Co radioactive source under the following conditions: 10% efficiency ORTEC GAMMA-X PLUS detector, 8.5-V amplitude for the 1.33-MeV gamma-ray on the unipolar output.

SPECTRUM SHIFT† (Fig. 2.1.) Peak position typically shifts ± 0.01 8% for counting rates up to 100,000 counts/s, and ± 0.05 % for counting rates up to 200,000 counts/s. Measured on the 1.33-MeV line under conditions specified for SPECTRUM BROADENING.

DIFFERENTIAL INPUT Differential nonlinearity ± 0.012 % from -9 V to +9 V. Maximum input ± 10 V (dc plus signal). Common mode rejection ratio >1000.

PULSE PILE-UP REJECTOR

Threshold Automatically set just above noise level on fast amplifier signal. Independent of slow amplifier BLR threshold.

Minimum Detectable Signal Limited by detector and preamplifier noise characteristics.

Pulse Pair Resolution Typically 500 ns. Measured using the ⁶⁰CO 1.33-MeV gamma-ray under the following conditions: 10% efficiency germanium detector, 4-V amplitude for the 1.33-MeV gamma-ray at the unipolar output, 50,000 counts/s.

†Results may not be reproducible if measured with a detector producing a large number of slow-risetime pulses or having quality inferior to the specified detector.

2.2. CONTROLS AND INDICATORS

FINE GAIN Front-panel, 10-turn precision potentiometer with locking, graduated dial provides continuously variable, direct reading, gain factor from 0.5 to 1.5.

COARSE GAIN Front-panel, eight-position switch selects gain factors of 5, 10, 20, 100, 200, 500, and 1000.

SHAPING TIME Six-position switch on the front panel selects shaping times of 0.5, 1, 2, 3, 6, and 10 μ s for the pulse-shaping filter network.

UNI SHAPING Two-position locking toggle switch on the front panel selects either GAUSSIAN or TRIANGLE pulse shaping for the UNIPOLAR output.

INPUT Front-panel, four-position switch accommodates either + or - input polarities, and selects the differential (DIFF) or normal (NORM) input modes. In the NORM mode only the NORM input connector is used. In the DIFF mode the preamplifier signal cable is connected to the NORM input, and a cable having its center conductor connected to the preamplifier ground through an impedance matching resistor is connected to the DIFF REF input. The impedance matching resistor must match the output impedance of the preamplifier.

DIFFERENTIAL INPUT GAIN BALANCE A 20-turn potentiometer mounted on the PC board inside the module allows the gains of NORM and DIFF REF inputs to be matched for maximum common mode noise rejection in DIFF mode.

PZ AUTO/MAN SWITCH Locking toggle switch selects either the AUTO (automatic) or MAN (manual) pole-zero cancellation adjustment mode. Both modes permit PZ cancellation for preamplifier exponential decay time constants from 40 μ s to ∞ .

AUTO PZ BUTTON With PZ switch in AUTO PZ position, momentarily pressing AUTO PZ button turns on the BUSY LED and initiates automatic adjustment of the PZ cancellation circuit. BUSY LED turns off when adjustment is complete. Once completed, the PZ adjustment is held until the button is pushed again. Memory of the last PZ adjustment is protected against unforeseen power outage.

MANUAL PZ ADJUSTMENT 20-turn potentiometer on the front panel permits screwdriver adjustment of the PZ cancellation. The screwdriver-adjusted value is effective whenever the PZ switch is in the MAN (manual) position. For transistor reset preamplifiers or pulsed optical feedback preamplifiers, use manual PZ adjustment set fully counterclockwise.

LIMIT PUSHBUTTON Inserts a diode limiter in series with the front-panel UNIPOLAR output connector and test point. Prevents overload distortions in the oscilloscope when observing accuracy of the PZ adjustment on the more sensitive oscilloscope ranges.

BLR RATE A front-panel, three-position, locking, toggle switch selects the baseline restorer rate. PZ position offers lowest fixed rate, for adjusting PZ cancellation. AUTO position matches the rate of the PZ position at low counting rates, but increases the restoration rate as the counting rate rises. HIGH rate position is provided for suppressing low frequency interference.

PUR ACCEPT/REJECT LED Multicolor LED indicates percentage of pulses rejected because of pulse pile-up. LED appears green for 0-40%, yellow for 40-70%, and red for >70% rejection.

2.3. INPUTS

NORM Front-panel, BNC connector accepts preamplifier signals of either polarity with risetimes less than the selected SHAPING TIME, and exponential decay time constants from 40 μ s to ∞ . For the — INPUT switch setting, the input impedance is 1000 Ω on a coarse gain of 5, and 465 Ω at coarse gain settings ≥ 10 . For the + INPUT switch setting, the input impedance is 2000 Ω for a coarse gain of 5, and 1460 Ω for coarse gains ≥ 10 . Input is dc-coupled, and protected to ± 25 V.

LINEAR Rear-panel connector. Identical to NORM input.

DIFF REF Front-panel BNC connector is used for the preamplifier ground reference connection when operating in the differential input mode. Operative only with the INPUT switch in the DIFF mode. For the + DIFF INPUT switch setting, the input impedance is 1000 Ω on a coarse gain of 5, and 465 Ω at coarse gain settings ≥ 10 . For the — DIFF INPUT switch setting, the input impedance is 2000 Ω for a coarse gain of 5, and 1460 Ω for coarse gains ≥ 10 . Input dc-coupled; protected to ± 25 V.

INHIBIT Rear-panel BNC input connector accepts reset signals from transistor reset preamplifiers or pulsed optical feedback preamplifiers. Positive NIM standard logic pulses or TTL levels can be used. Logic is selectable as active high or active low via printed circuit board jumpers. INHIBIT input initiates the protection against distortions caused by the preamplifier reset. This includes turning off the baseline restorers, monitoring the negative overload recovery at the unipolar output, and generating PUR (reject) and BUSY signals for the duration of the overload. The PUR and BUSY logic pulses are used to prevent analysis and correct for the reset deadtime in the associated ADC or multichannel analyzer.

2.4. OUTPUTS

UNIPOLAR, UNI Front- and rear-panel BNC connectors provide positive, unipolar, shaped pulses with a linear output range of 0 to +10 V. Front-panel output impedance < 1 Ω . Rear-panel output impedance selectable for either < 1 Ω or 93 Ω using a printed circuit board jumper. Outputs are dc-restored to 0 ± 5 mV and short-circuit protected.

BIPOLAR, Bi Front- and rear-panel BNC connectors provide bipolar shaped pulses with the positive lobe leading. The linear output range is 0 to ± 10 V. Front- panel output impedance $< 1 \Omega$. Rear-panel output impedance selectable for either $< 1 \Omega$ or 93Ω using a printed circuit board jumper. Baseline between pulses has a dc level of 0 ± 10 mV. Short-circuit protected.

CRM The Count Rate Meter output has a rear-panel BNC connector and provides a 250-ns-wide, +5-V logic signal for every linear input pulse that exceeds the pile-up inspector threshold. Output impedance is 50Ω .

BUSY Rear-panel BNC connector provides a +5-V logic pulse for the duration that the linear signals exceed the positive or negative baseline restorer thresholds, or the pile-up inspector threshold, or for the duration of the INHIBIT input signal. Useful for deadtime corrections with an associated ADC or multichannel analyzer. Positive NIM standard logic pulse is selectable as active high or active low via a printed circuit board jumper. Output impedance is 50Ω .

PUR Pile-Up Reject output is a rear-panel, BNC connector. Provides a +5-V NIM standard logic pulse when pulse pile-up is detected. Output also present for a pulsed reset preamplifier during reset, and reset overload recovery. Output pulse is selectable as active high or active low by means of a printed circuit board jumper. Output impedance is 50Ω . Used with an associated ADC or multichannel analyzer to prevent analysis of distorted pulses (Fig. 2.2).

PREAMP Rear-panel standard ORTEC connector (Amphenol 17-10090) provides power for the associated preamplifier. Mates with power cords on all standard ORTEC preamplifiers.

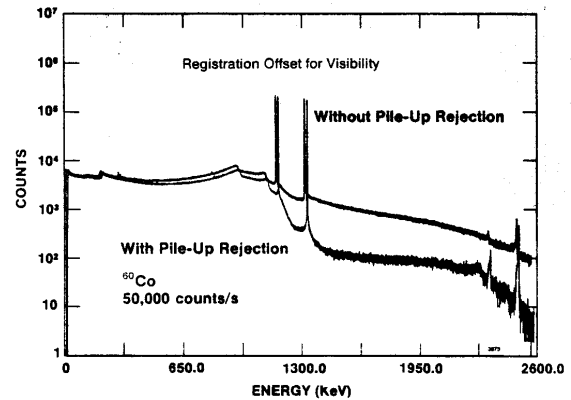


Fig. 2.2. Demonstration of the Effectiveness of the Pile-Up Rejector in Suppressing the Pile-Up Spectrum.

2.5. ELECTRICAL AND MECHANICAL

POWER REQUIRED The Model 672 derives its power from a NIM Bin supplying ± 24 V and ± 12 V, such as the ORTEC Model 4001A/4002A Bin/Power Supply. The power required is +24 V at 90 mA, -24 V at 170 mA, +12 V at 330 mA, and -12 V at 190 mA.

WEIGHT

Net 2.3 kg (5.0 lb).

Shipping 3.6 kg (8.0 lb).

DIMENSIONS NIM-standard double-width module 6.90 X 22.13 cm (2.70 X 8.714 in.) per DOE/ER-0457T.

3. INSTALLATION

3.1. POWER CONNECTION

The 672 operates on power that must be provided by a NIM-standard bin and power supply such as the ORTEC 4001/4002 series. Convenient test points on the power supply control panel should be used to check that the dc voltage levels are not overloaded. The bin and power supply is designed for relay rack mounting. If the equipment is rack mounted, be sure that there is adequate ventilation to prevent any localized heating of the components that are used in the 672. The temperature of the equipment mounted in racks can easily exceed the maximum limit of 50°C unless precautions are taken.

3.2. PREAMPLIFIER CONNECTION

The Preamp connector of this amplifier is directly compatible with ORTEC preamplifiers as well as with standard Aptec, Canberra, PGT, and Tennelec (serial numbers greater than 2000) preamplifiers. Preamplifier power at +24 V, -24 V, +12 V, and -12 V is available through the Preamp connector on the rear panel.

When a BNC cable longer than ten feet is used to connect the preamplifier output to the amplifier input, the characteristic impedance of the cable should match the impedance of the preamplifier output. All ORTEC preamplifiers contain series terminations that are either 93 Ω or variable; coaxial cable type RG-62/U or RG-71/U is recommended.

3.3. PULSED RESET PREAMPLIFIERS AND INHIBIT IN CONNECTION

The 672 Amplifier is directly compatible with most pulsed reset preamplifiers such as the ORTEC TRP (Transistor Reset Preamplifier) Series. The amplifier automatically senses preamplifier resets and gates off the amplifier's baseline restorer. Preamplifier inhibit signals are not required for proper amplifier operation; however, since the

preamplifier resetting process is nonlinear by nature, spurious phantom peaks may show up in the spectra if the inhibit signal from the preamplifier is not used.

INHIBIT IN CONNECTION Connection of the PREAMPLIFIER INHIBIT OUT signal to the rear-panel INHIBIT IN connector will result in the system being disabled during the reset period and thus avoid spurious peaks in the spectra. Preamplifiers with an Inhibit time switch such as ORTEC PLUS Detector with series 132 Preamplifier can be set to position "1", which is the shortest preamp inhibit blocking time.

PZ SETTING The Amplifier's Manual PZ control should be set fully counterclockwise (CCW) when used with a pulsed reset preamplifier. The front-panel control is marked with an arrow labeled TRP to indicate the direction to turn the control. The Auto PZ should not be used with Pulsed Reset Preamplifiers.

3.4. CONNECTION OF TEST PULSE GENERATOR THROUGH A PREAMPLIFIER

THROUGH A PREAMPLIFIER The satisfactory connection of a test pulse generator such as the ORTEC 419 or 448 Pulse Generator or equivalent depends primarily on two considerations: the preamplifier must be properly connected to the 672 as discussed in Sections 3.2 and 3.3, and the proper signal simulation must be applied to the preamplifier. To ensure proper input signal simulation, refer to the instruction manual for the particular preamplifier being used.

DIRECTLY INTO THE 672 The ORTEC test pulse generators are designed for direct connection. When any one of these units is used, it should be terminated with a 100- Ω terminator at amplifier input or be used with at least one of the output attenuators set at In.

SPECIAL CONSIDERATIONS FOR POLE-ZERO CANCELLATION When a tail pulser is connected directly to the amplifier input, the Pole-Zero should be adjusted. See Section 4.3 for the pole-zero adjustment. If a preamplifier is used and a tail pulser is connected to the preamplifier test input, it is not possible to adjust the pole-zero for both the preamplifier pole and the pole from the pulser tail.

3.5. SHAPING CONSIDERATIONS

The Shaping Time switch on the front panel of the 672 can be set to select time constants in steps of 0.5, 1, 2, 3, 6, and 10 μs . Choice of triangular and Gaussian filters doubles the time constants available for optimum resolution. Triangular shaping will usually give better results. The choice of the proper shaping time is generally a compromise between operating at a shorter time constant for accommodation of high counting rates and operating with a longer time constant for a better signal-to-noise ratio. Since the full amplitude of the preamplifier output pulse must be preserved, the peaking time (measurement time) must be large compared to preamplifier output pulse risetime. The amplifier shaping time should be greater than five times the charge collection time of the detector. Use the detector manufacturer's suggested shaping times as a starting point and adjust the shaping as your needs for resolution versus count rate vary.

GERMANIUM DETECTORS Shaping times for high-purity germanium (HPGe) detectors will vary from 1 to 6 μs using the Unipolar output, depending on the size, configuration, and charge collection time of the specific detector and preamplifier. Coaxial detectors have significant variations in charge collection times due to their large volumes. Compromises must often be made since the shaping time that will give the best resolution will usually be longer than the optimum time needed for the best throughput at high counting rates.

Planar detectors require shaping times in the range of 3 to 10 μs for optimum resolution. Lithium-drifted silicon detectors, Si(Li), have similar shaping time requirements.

SILICON CHARGED PARTICLE DETECTORS

These detectors have very fast risetimes on the order of 10 ns or less. A Unipolar output and a 0.5- to 2- μs shaping time will generally provide optimum resolution.

SCINTILLATION DETECTORS The energy resolution of scintillation counters depends largely on the scintillator and photomultiplier, and therefore a shaping time of five times the decay-time constant of the scintillator is a reasonable choice. For NaI detectors that have a decay time constant of about 230 ns, the optimum shaping time is 1 μs . The Bipolar output can be used to reduce overload effects and microphonics without sacrificing resolution.

GAS PROPORTIONAL COUNTERS Proportional counters have both short and long components in their charge collection times. The components typically fall in the 0.5- to 5- μs range, and lead to variable amounts of preamplifier output signal being lost as the amplifier shaping time constant is changed. Selection of longer shaping times ($>2 \mu\text{s}$) helps to minimize the problem caused by long risetimes. Due to the multiple components in the charge collection time, the correct pole-zero cancellation is not possible. This will often cause an undershoot if the Unipolar output issued. Bipolar shaping can be used to reduce this effect with little change in the resolution.

3.6. LINEAR OUTPUT CONNECTIONS

Since the 672 Unipolar output is normally used for spectroscopy, the 672 is designed with a great amount of flexibility for the pulse to be interfaced with an analyzer. To minimize spectrum distortion at medium and high counting rates, the unipolar output incorporates a high-performance, gated baseline restorer with automatic setup. Automatic positive and negative noise discriminators ensure that the baseline restorer operates only on the true baseline between pulses in spite of changes in the noise level. For pulse-height analysis the Unipolar output must be directly connected to the input of a multichannel analyzer.

The bipolar output, with its symmetry about the baseline, can be used for cross-over timing or may be preferred for spectroscopy when operating into ac-coupled systems at high counting rates. Typical system block diagrams for a variety of experiments are described in Section 4.

3.7. PILE-UP REJECTION USING PUR OUTPUT

The PUR (Pile-Up Reject) output on the rear panel is used at the gate or pile-up reject input of a multichannel analyzer to suppress pile-up in the recorded spectrum. The fast amplifier in the pile-up rejector includes a gated baseline restorer with an automatic noise discriminator to eliminate the need for any operator adjustments. When pileup occurs, a logic true pulse is generated which lasts until the unipolar output returns to the baseline, normally a width of six times the shaping time. If used with a pulsed reset preamplifier, this output also includes a reject during the reset and recovery interval.

3.8. LIVETIME CORRECTION USING BUSY OUTPUT

The signal from the rear-panel Busy output connector provides a nominally +5 V logic pulse for the duration that the Unipolar output pulse exceeds the baseline restorer threshold or pile-up inspector threshold or when the external INHIBIT IN is true. For livetime correction, Busy should be connected to the Busy In connector on the MCA. For optimal livetime correction with ORTEC analyzers like the ADCAM[®], an internal jumper in the amplifier should be set to match the unipolar, triangular, or Gaussian mode. The output is internally jumper selectable as active low or active high. It is shipped as active high.

3.9. INPUT COUNT RATE USING CRM OUTPUT

A positive logic pulse is generated for each 672 input pulse that exceeds the pile-up inspector threshold level. The pulses are available through the CRM (Count Rate Meter) output on the rear panel and are intended for use in a count rate meter or counter to monitor the true input count rate into the amplifier.

4. OPERATING INSTRUCTIONS

4.1. INITIAL TESTING AND OBSERVATION OF PULSE WAVEFORMS

Refer to Section 6 for information on testing performance and observing waveforms using a pulser. Figure 4.1 shows some typical unipolar Gaussian, unipolar triangular, and bipolar output waveforms.

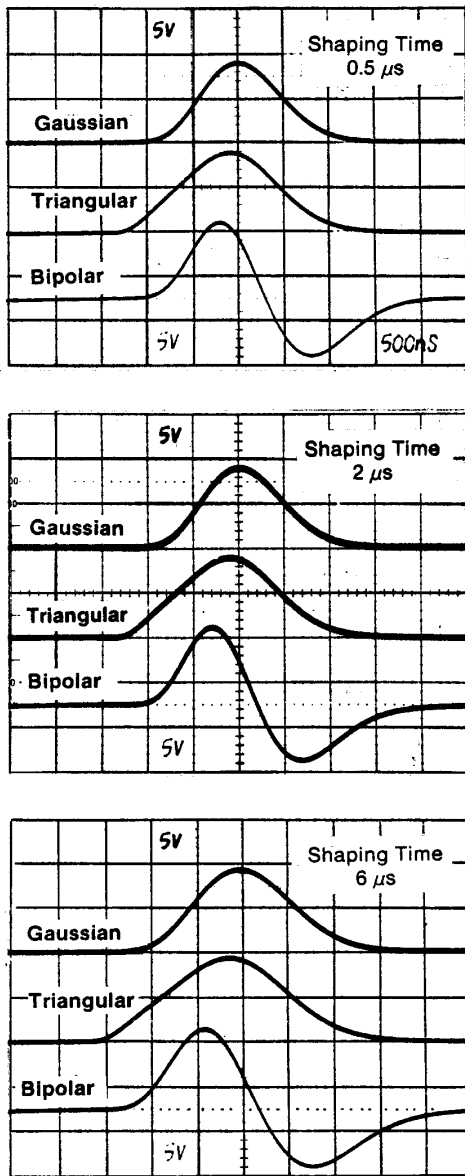


Fig. 4.1. Typical Effects of Shaping-Time Selection on Gaussian, Triangular, and Bipolar Output Waveforms.

4.2. STANDARD SETUP PROCEDURES

a. Connect the detector, preamplifier, high-voltage power supply, and amplifier into a basic system and connect the amplifier unipolar output to an oscilloscope. Connect the preamplifier power cable to the Preamp power connector on the rear panel of the 672. Turn on power in the bin and power supply and allow the electronics of the system to warm up and stabilize.

A block diagram of a typical ORTEC gamma-ray spectroscopy system is shown in Figure 4.2.

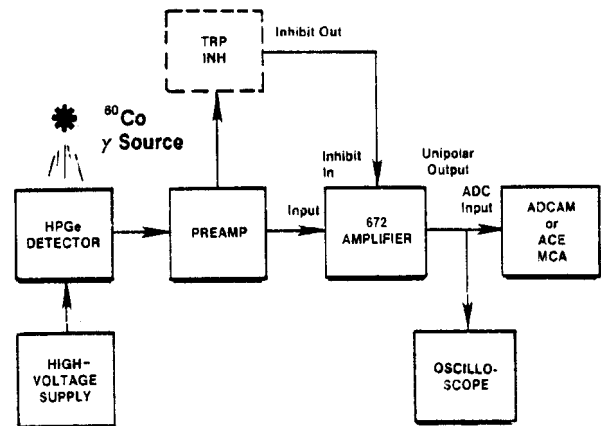


Fig. 4.2. Typical Gamma-Ray Spectroscopy System.

b. Set the 672 controls initially as follows:

Shaping	3 or 6 μ s
Unipolar Mode	Triangle
Coarse Gain	20
Fine Gain	1.00
BLR Rate	PZ
Polarity	Norm, Match preamplifier output polarity
PZ	Auto

c. Use a ^{60}Co calibration source; set about 25 cm from the active face of the detector. The unipolar output pulse from the 672 should be about 8 V, using a detector that has a preamp with a conversion gain of 300 mV/MeV.

d. Readjust the Gain control so that the higher peak from the ^{60}Co source (1.33 MeV) provides an amplifier output at about 9 V.

4.3. POLE-ZERO ADJUSTMENTS FOR RESISTIVE-FEEDBACK PREAMPLIFIER

The pole-zero adjustment is critical for good performance at high count rates in unipolar operation and for correct operation of the BLR circuit. This adjustment should be checked carefully for the best possible results. Whenever the shaping time is changed, the pole-zero must be adjusted. The bipolar output resolution is not as sensitive to misadjusted PZ, but it is important for recovery from very large overload pulses. When using a transistor reset-type preamplifier, the Manual PZ Mode should be used and should be set to full counterclockwise.

AUTO POLE-ZERO MODE This mode should give good results for most detectors for low to medium count rates. For very high count rates when the count rate times the dead time per pulse (T_d) is greater than 0.8, resolution may sometimes be improved with some detector-preamp combinations by manually adjusting the PZ. The dead time, T_d , for the 672 is about nine times the front-panel shaping time. Memory of the last PZ adjustment is protected against unforeseen power outages.

a. Using the system detector and a radioactive source as an input, adjust the source distance for an input count rate of between 1 and 10 kHz.

b. Select Auto PZ, depress the push button and wait a few seconds. The red LED will go on and then turn off when the unipolar output pulse has been pole-zero compensated.

No oscilloscope is required for this procedure. However, the quality of the PZ can be observed with an oscilloscope by increasing the scope input

sensitivity to 10 or 20 mV/cm with the LIMIT button pushed. The LIMIT push-button switch inserts a diode limiter in series with the front-panel Unipolar Output connector to prevent overloading the input of the oscilloscope.

NOTE: The Auto PZ must be reset after power down.

MANUAL POLE-ZERO MODE

a. Adjust the radiation source spacing from the detector to provide a count rate between 1 and 10 kHz.

b. Observe the unipolar output with an oscilloscope. Increase the scope input sensitivity to 20-100 mV per vertical division. Depress the front-panel LIMIT push-button to limit the voltage applied to the oscilloscope. Adjust the PZ adjust control so that the trailing edge of the pulses returns to the baseline without overshoot or undershoot (Fig. 4.3). A slight bias toward an undershoot often gives the best results.

The oscilloscope used must be dc-coupled and must not contribute distortion in the observed waveforms. Oscilloscopes such as Tektronix models 465, 475, and 7904 will overload for a 10-V signal when the vertical sensitivity is <100 mV/Div. The LIMIT push-button switch inserts a diode limiter in series with the front-panel Unipolar Output connector to prevent overloading the input of the oscilloscope.

USING SQUARE WAVE THROUGH PREAMPLIFIER TEST INPUT

A more precise manual pole-zero adjustment of the amplifier can be obtained by using a squarewave signal as the input to the preamplifier. Many oscilloscopes include a calibration output on the front panel, and this is a good source of square wave signals at a frequency of about 1 kHz. The amplifier differentiates the signal from the preamplifier so that it generates output signals of alternate polarities on the leading and trailing edges of the square wave input signal, and these can be compared as shown in Fig. 4.4 to achieve excellent pole-zero cancellation.

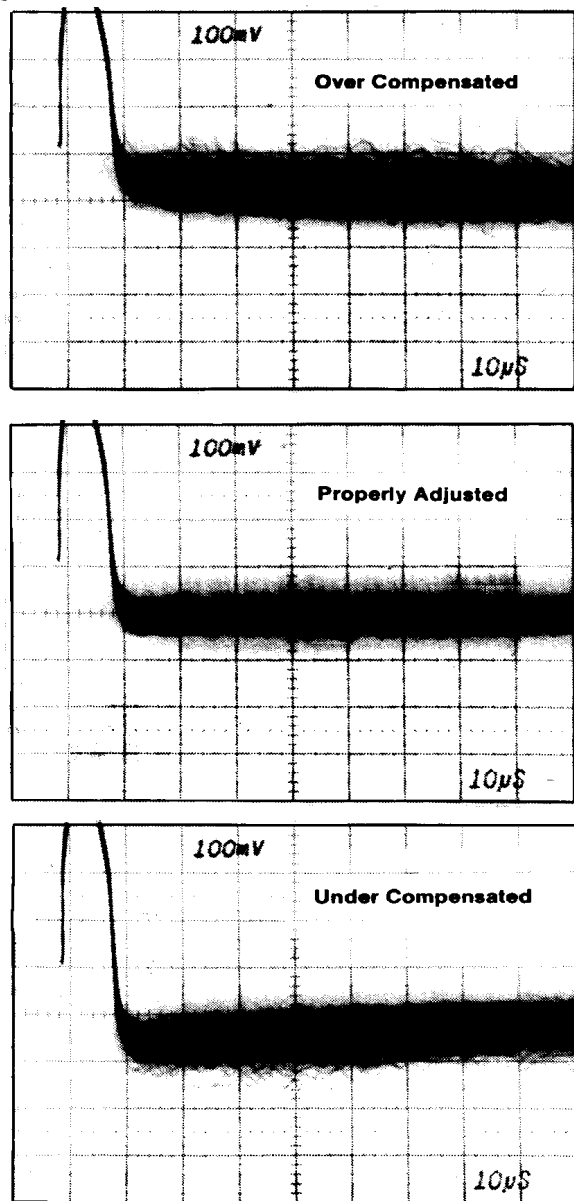


Fig. 4.3. Typical Waveforms Illustrating Pole-Zero Adjustment Effects; Oscilloscope Trigger, Busy Output; ^{60}Co Source with 1.33-MeV Peak Adjusted ~ 9 V; Count Rate, 3 kHz; Shaping Time Constant, 2 μs .

Use the following procedure:

- Remove all radioactive sources from the vicinity of the detector. Set up the system as for normal operation, including detector bias.
- Set the amplifier controls as for normal operations; this includes gain, shaping, and input polarity.

c. Connect the source of 1 kHz square waves through an attenuator to the Test input of the preamplifier. Adjust the attenuator so that the amplifier output amplitude is 8 to 10 volts.

d. Observe the unipolar output of the amplifier with an oscilloscope triggered from the amplifier Busy output. Adjust the PZ control for proper response according to Fig. 4.4. Depress the LIMIT push-button on the 672 while observing the adjustment on the oscilloscope display.

Figure 4.4A shows the amplifier output as a series of alternate positive and negative shaped pulses. In B, C, and D of this figure the oscilloscope was triggered to show both positive and negative pulses simultaneously. These pictures show more detail to aid in proper adjustment.

4.4. BASELINE RESTORER (BLR) SETTING

To minimize spectrum distortion at medium and high counting rates, the unipolar output incorporates a high-performance, gated, baseline restorer with several levels of automation. Automatic positive and negative noise discriminators ensure that the baseline restorer operates only on the true baseline between pulses in spite of changes in the noise level. No operator adjustment of the baseline restorer is needed when changes are made in the gain, the shaping time constant, or the detector characteristics. Negative overload recovery from the reset pulses generated by transistor reset preamplifiers and pulsed optical feedback preamplifiers is also handled automatically to eliminate the need for operator adjustments. A monitor circuit gates off the baseline restorer and provides a reject signal for a multichannel analyzer until the baseline has safely recovered from the overload.

BLR RATE For making either a manual or automatic pole-zero adjustment, the PZ position is selected. This position can also be used where the slowest baseline restorer rate is desired.

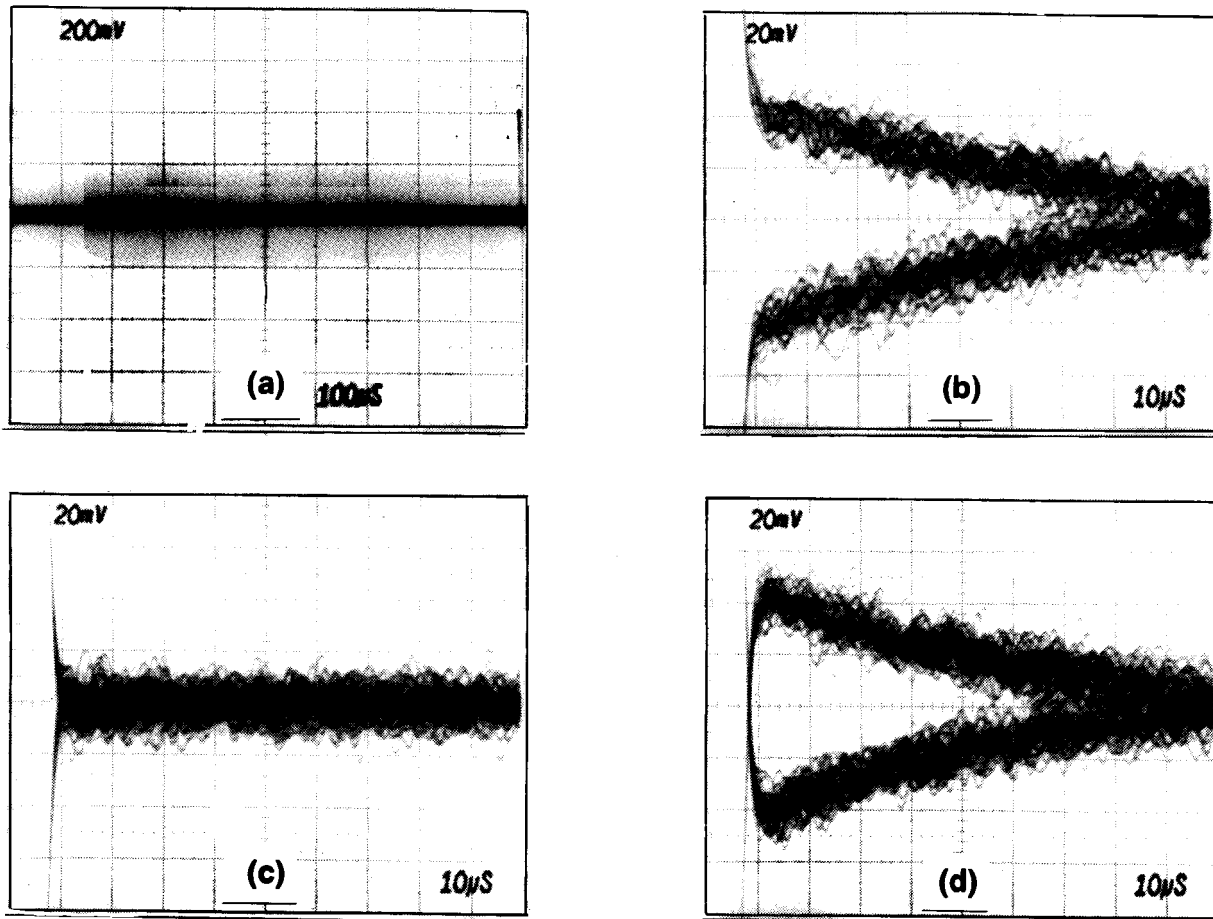


Fig. 4.4. Pole-Zero Adjustment Using a Square Wave Input to the Pre-amplifier. (a) PZ properly adjusted; slow trigger to separate pulses. (b) Overcompensated; fast trigger to superimpose pulses. (c) Properly adjusted; pulses superimposed. (d) Undercompensated; pulses superimposed.

With the BLR Rate set to AUTO, the BLR is automatically set for optimum performance throughout the usable input range for the shaping selected.

The HIGH rate can be used for situations where low or medium frequency noise interference is present and is independent of the counting rate. The HIGH rate setting is normally not used since there will be a small loss of resolution due to increased noise when used in high resolution systems.

4.5. INTERNAL CONTROLS

These controls are on the printed wiring board (PWB) and can be accessed by removing the right

side cover. Figure 4.5 shows the location of these controls.

BAL (DIFFERENTIAL INPUT GAIN BALANCE)

Internal PWB 20-turn screwdriver potentiometer allows maximization of noise rejection when using differential input. See Section 4.6.

UNI-OUT (UNIPOLAR Z_{OUT}) Jumper plug, W1, provides Z_{OUT} , $\leq 1 \Omega$ or $\sim 93 \Omega$) for the rear-panel Unipolar output. Shipped in the $93\text{-}\Omega$ position.

BI-OUT (BIPOLAR Z_{OUT}) Jumper plug, W2, provides Z_{OUT} , $\leq 1 \Omega$ or $\sim 93 \Omega$) for the rear-panel Bipolar output. Shipped in the $93\text{-}\Omega$ position.

BUSY/BUSY Jumper plug, W3, allows the Busy output to be a positive true or negative true logic signal. Shipped in BUSY (positive true) position.

PUR/PUR Jumper plug, W5, allows the Pile-Up Reject (PUR) output to be a positive true or negative true logic signal. Shipped in PUR (positive true) position.

INH/INH Jumper plug, W6, allows the inhibit in input to accept either positive true or negative true logic signals. Shipped in INH (positive true) position.

TRI/GAUSS Jumper plug, W7, allows optimal livetime correction when used with ORTEC analyzers like the ADCAM® by connecting the BUSY output to the analyzer Busy In as described in Section 3.8. The jumper should be set to match the Unipolar Mode, TRI for Triangle and GAUSS for Gaussian. Shipped in TRI position.

4.6. DIFFERENTIAL INPUT MODE

When long connecting cables are used between the detector and preamplifier input, noise induced in the

cable by the environment can be a problem. The differential input mode can be used with paired cables from the preamplifier to suppress the induced noise.

BAL (DIFFERENTIAL INPUT GAIN BALANCE)

The BAL potentiometer is used to adjust the gain balance between the positive and negative inputs and to adjust the balance between the NORM and DIFF REF inputs when the differential input mode is used. The initial adjustment of Gain Balance is made by providing the same input to both the NORM and DIFF REF inputs. This can be accomplished by using a BNC "T" connector to feed the input signal on the NORM input to the DIFF REF input. Set the amplifier gain to maximum. Connect an oscilloscope to the unipolar output. While observing the signal on the oscilloscope, used a small screwdriver to adjust the Gain Balance (internal adjustment has been factory set, Fig. 4.5) potentiometer until the display on the oscilloscope shows minimum signal. Remove the BNC "T" connector when the adjustment is complete, and the positive and negative gains will be matched for use with NORM input.

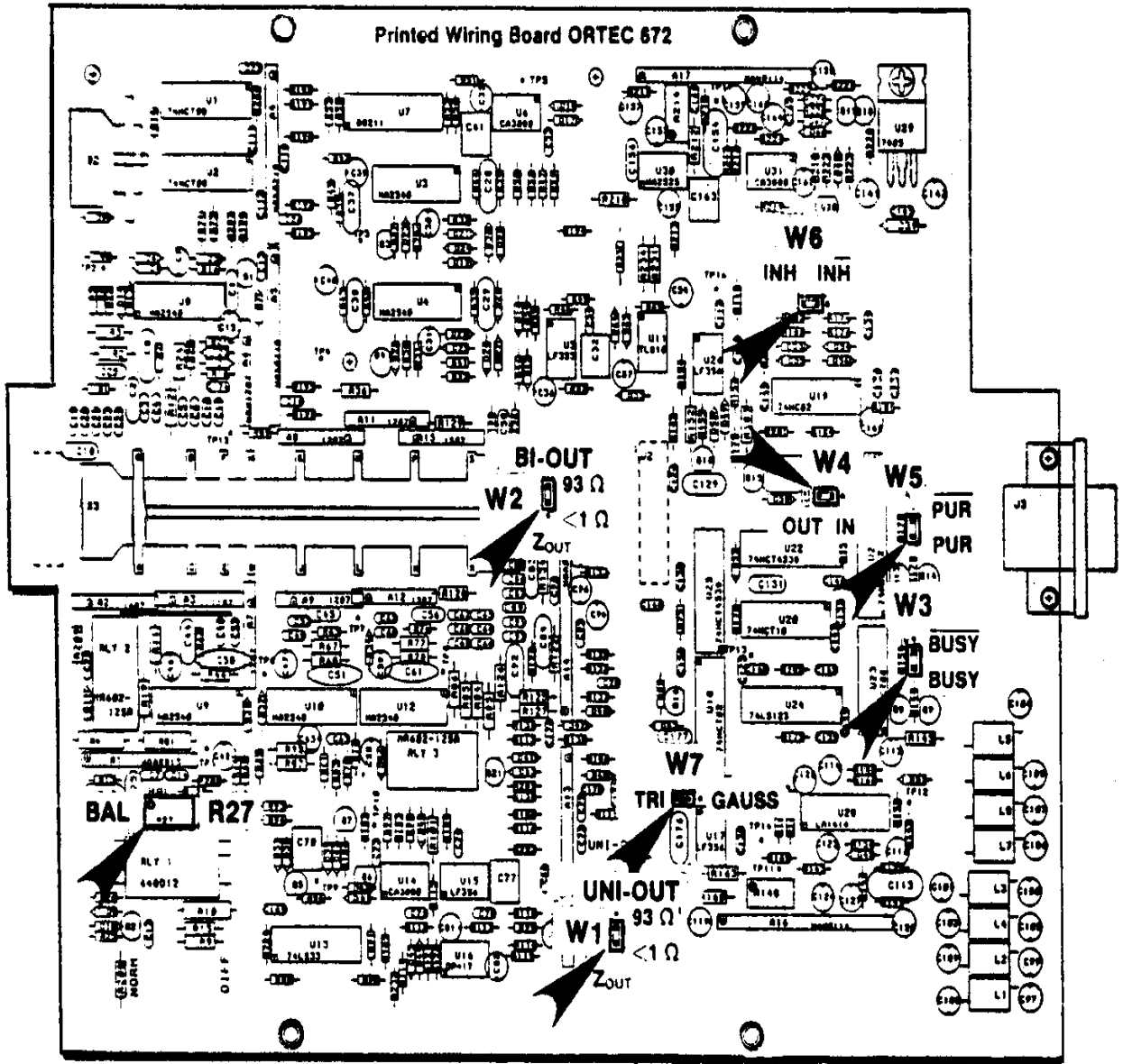


Fig. 4.5. Position of Internal Controls.

If the differential input mode is being used, connect the differential input cable to the DIFF BNC connector on the front panel. Adjust BAL potentiometer until there is minimum noise around the baseline of the output signal. If there is a problem in getting minimum noise, repeat the initial procedure with the BNC "T" and the adjustment.

DIFF INPUT SIGNAL The DIFF input signal or phantom is used only in the differential input mode. The normal preamp output is connected to the NORM input with the amplifier input polarity set to match this signal. A second output cable must be added to the preamplifier with its center, signal pin connected to the preamplifier ground with the same value as the normal preamp output series resistor (usually 93.1 or 51 Ω).

Many ORTEC preamplifiers have two Energy outputs, each with a 93.1- Ω series resistor. For differential operation, one output is connected to the amplifier NORM input. The second output is modified by connecting the preamplifier end of the series 93.1- Ω resistor to ground within the preamp (soldering may be necessary). This second output should be properly marked and connected to the DIFF REF input. Both cables should be the same length and be run next to each other.

4.7. SYSTEM THROUGHPUT

To achieve the desired results in high-rate energy spectroscopy, the experimenter must consider not only the input rate, but also the unpiled-up output rate. The unpiled-up output rate is determined by the processing time of the shaping amplifier, the pile-up inspection time, and the input rate. For semi-Gaussian time-invariant filter amplifiers, the unpiled-up output rate is theoretically given by²

$$r_o = r_i \exp(-T_D r_i) \quad (1)$$

where r_o is the unpiled-up output count rate, r_i is the input count rate, and T_D is the deadtime or effective processing time of the amplifier. The value of T_D is equal to the sum of the effective amplifier pulse

width, T_w , and the time-to-peak of the amplifier output pulse, T_p . The type of deadtime in the shaping amplifier is referred to as extending deadtime since a second event arriving before the end of the initial deadtime extends the deadtime by an additional amplifier output pulse width, T_w , from the occurrence of the second pulse.

A normalized plot of Equation (1) is shown as the solid line in Fig. 4.6. The maximum mean output rate equals $1/T_D \exp(1)$ and occurs when the mean input rate equals $1/T_D$. At this maximum output rate the deadtime losses are 63.2%. For input count rates exceeding $1/T_D$ the unpiled-up output rate decreases. When using a pile-up inspection circuit, the value of T_D is given either by the sum of T_w and T_p , or by the sum of T_p and the pile-up inspection time, whichever is larger.

Spectroscopy systems also have a deadtime that is caused by the digitizing time of the Analog-to-Digital Converter (ADC). This deadtime is a non-extending deadtime since events arriving during the digitizing time are ignored. For non-extending deadtime the output rate is given by²

$$r_o = \frac{r_i}{1 + r_i T_D} \quad (2)$$

where T_D is the digitizing time for the ADC and is designated T_M in Equation (3). This relationship is shown as the dashed line in Fig. 4.6. The maximum obtainable output count rate is $1/T_D$ and occurs at $r_i = \infty$.

When the ADC is considered as part of the spectroscopy system, the deadtimes of the amplifier and ADC are in series. The combination of the extending deadtime of the amplifier followed by the non-extending deadtime of the ADC is given by²

$$r_o = \frac{r_i}{\exp[r_i(T_w + T_p)] + r_i[T_M - (T_w - T_p)] U[T_M - (T_w - T_p)]} \quad (3)$$

²R. Jenkins, R.L. Gould, and D.A. Gedcke, *Quantitative X-Ray Spectroscopy*, Marcel and Dekker, Inc., New York, (1980).

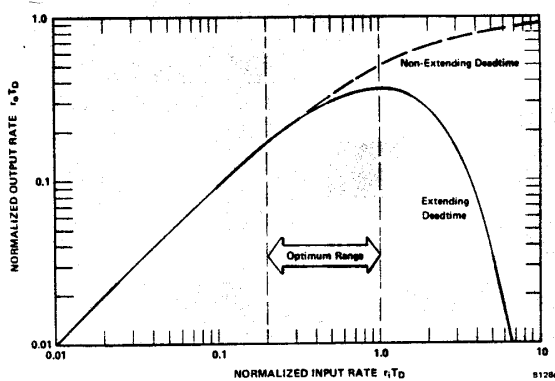


Fig. 4.6. Plot of Normalized Output Rate as a Function of Normalized Input Rate for Spectrometers with Simple Deadtime.

where $U[T_M - (T_W - T_P)]$ is a unit step function that changes value from 0 to 1 when T_M is greater than $(T_W - T_P)$. Equation (3) reduces to Equation (1) when T_M is less than $(T_W - T_P)$.

A plot of the unpiled-up amplifier output rate as a function of input rate for six values of shaping time is shown in Fig. 4.7. The measured deadtime, T_D , is shown for each shaping time constant. The maximum value of the unpiled-up output rate increases with decreasing values of shaping time constant. A set of throughput curves will remain nearly unchanged for a given amplifier for various energy ranges, detector types, and sizes.

The advantage of shorter shaping time constants to achieve higher output count rates is clearly shown in Fig. 4.7. However, shorter time constants also result in increased noise and increased charge collection time effects. Under worst case conditions, the noise increases inversely as the square root of the ratio of shaping time constants. The increase in the total energy resolution is the noise contribution combined in quadrature with the statistical contribution of the detector at the energy of interest. Consequently, the percentage of degradation in energy resolution can be much less than the percentage increase in noise.

4.8. CHARGE COLLECTION OR BALLISTIC DEFICIT EFFECTS

Charge collection distances in large-volume HPGe detectors are often 3 cm or more, resulting in

charge collection times exceeding 300 ns.^{3,4,5} These charge collection times are due to the transit time of the holes and the electrons in germanium and are not due to defects in the detector. Fig. 4.8(a) shows some typical current pulse waveforms from a 140-cm³ 28% efficient HPGe detector. These current pulse waveforms were obtained using the simple differentiation circuit shown in Fig. 4.8(b), which has a 15-ns time constant. The current pulses range in duration from 100 ns to greater than 350 ns. Pulses having equivalent total charge but different durations produce different output pulse heights when processed by a charge-sensitive preamplifier and a semi-Gaussian filter amplifier. This results in the distortion of the spectrum in direct proportion to the pulse amplitude or energy. This distortion is most pronounced at short shaping time constants. Figure 4.9(a) shows a portion of a spectrum obtained with a 10% efficient HPGe detector at 2- μ s shaping time, using the 1.33-MeV line of ⁶⁰Co. An equivalent spectrum using a 0.5- μ s shaping time is shown in Fig. 4.9(b) and is significantly distorted.

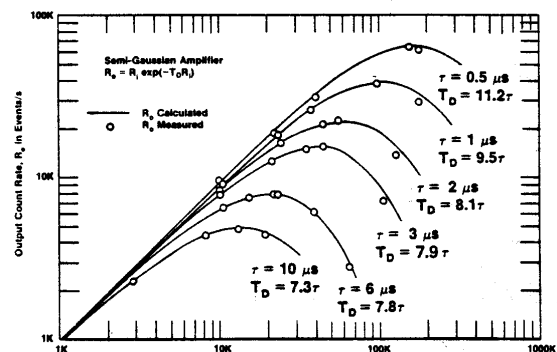


Fig. 4.7. Plot of the Unpiled-Up Amplifier Output Rate as a Function of Input Rate for Six Values of Shaping Time Constants.

³E. Sakai, "Charge Collection in Coaxial Ge(Li) Detectors," *IEEE Trans. Nucl. Sci.*, **NS-15**, 310, (1968).

⁴E. Sakai, T.A. McMath, and R.G. Franks, "Further Charge Collection Studies in Coaxial Ge(Li) Detectors," *IEEE Trans. Nucl. Sci.*, **NS-16**, 68, (1968).

⁵T.H. Becker, E.E. Gross, and R.C. Trammell, "Characteristics of High-Rate Energy Spectroscopy Systems with Time-invariant Filters," *IEEE Trans. Nucl. Sci.*, **NS-28**, 1, (1981).

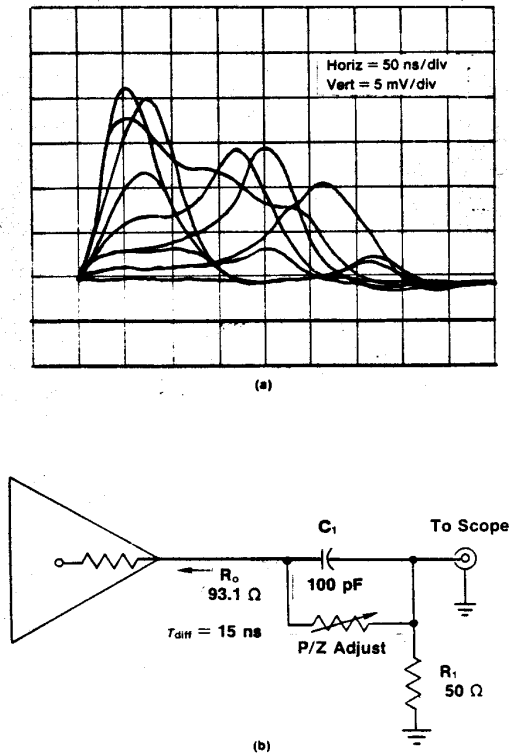


Fig. 4.8. Charge Collection Effect Waveforms. (a) Typical current Pulse Waveforms for a 28% Efficient HPGe Detector, and (b) the Simple Differentiation Circuit Used to Obtain the Current Waveforms.

Charge collection time effects are of significant importance when using large-volume Ge detectors at high energy. The performance of two HPGe detectors is compared in Fig. 4.10 at two different energies. When using the 122-keV line of ^{57}Co , the principal cause of resolution degradation with decreased shaping time constant is the increase in noise. However, when using the 1.33-MeV line of ^{60}Co , the significant degradation in resolution is due to charge collection effects. The calculated resolution for the 10% detector at 1.33 MeV is shown as the dashed line in Fig. 4.10, and indicates approximately 2.0 keV FWHM at a 0.5- μs shaping time constant. The measured resolution under these test conditions was 7.2 keV, indicating that charge collection effects dominate. In Fig. 4.10, charge collection effects begin to appear at time constants less than 3 μs .

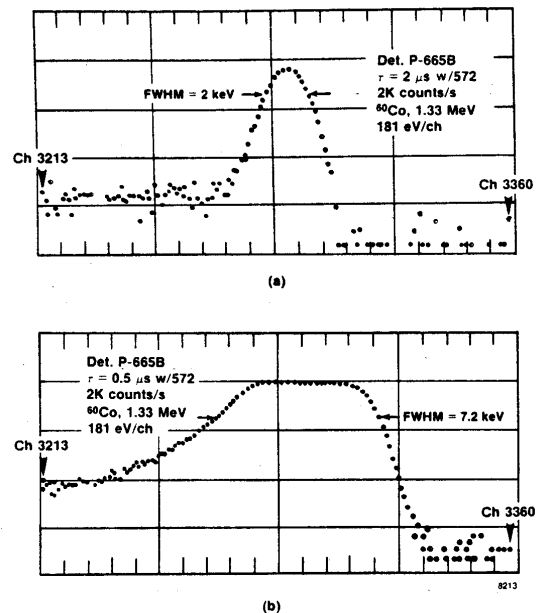


Fig. 4.9. Charge Collection Effect Spectrum. Logarithmic Display of Spectrum Taken with a 10% Efficient HPGe Detector for the 1.33 MeV ^{60}Co Line. (a) A 2- μs Shaping Time Constant and (b) a 0.5- μs Shaping Time Constant.

4.9. PILE-UP REJECTOR (PUR) AND LIVETIME CORRECTOR

An efficient pile-up rejector is incorporated in the amplifier to suppress the spectral distortion which is caused by pulses piling up on each other at high counting rates. High counting rate for pile-up is dependent on the dead time per pulse, T_D , and hence the selected shaping time. T_D is 9 times the front-panel shaping time, T_C . High count rate for the PUR is when the normalized count rate $R_i T_D > 0.5$, where R_i is the amplifier input rate (see Fig. 4.6). For example, for 6- μs shaping R_i is 9 kHz and for 2- μs shaping, R_i is 28 kHz. Amplifier throughput for this condition using Equation (1) in Section 4.7 is 60% of the input rate. A multicolor pile-up rejector LED is included on the front panel to indicate the throughput efficiency of the amplifier. At low counting rates (pulse pile-up losses <40%) the LED flashes with a green color. At moderate counting rates the color changes to yellow. The color changes to red at high counting rates when the pulse pile-up losses are >70%.

The fast amplifier in the pile-up rejector includes a gated baseline restorer with its own automatic noise discriminator to eliminate the need for any operator

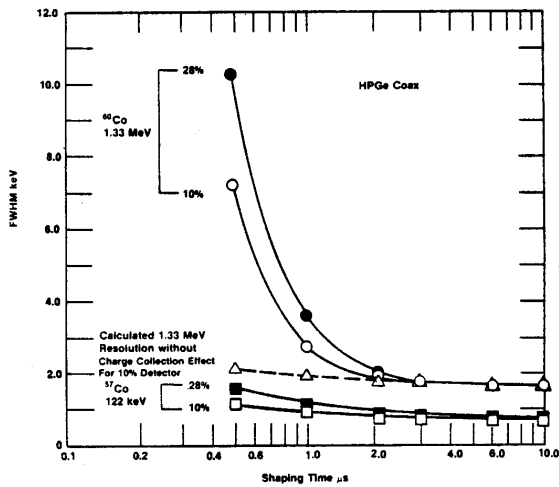


Fig.4.10. Energy Resolution FWHM as a Function of Amplifier Shaping Time Constant for a 10% HPGe Detector and a 28% HPGe Detector for the 122-keV ⁵⁷Co Line and the 1.33-MeV ⁶⁰Co Line.

adjustments. This function is also protected against negative overloads from pulsed reset preamplifiers. The PUR (pile-up reject) output logic pulse can be used at the gate or reject input of a multichannel

analyzer to suppress pile-up in the recorded spectrum.

The block diagram for a gamma-ray spectroscopy system with pile-up rejection and live time correction is shown in Fig. 4.11.

FOR A RESISTIVE FEEDBACK PREAMP, CONNECT:

a. Inhibit pulse from PUR to ADC PUR or ADC anticoincidence input.

b. Livetime correction signal (Busy output) to the ADC Busy In.

ADDITIONAL CONNECTION FOR TRP (Transistor Reset Preamplifiers) Shown in dotted lines.

c. Inhibit Output from TRP to the amplifier inhibit In.

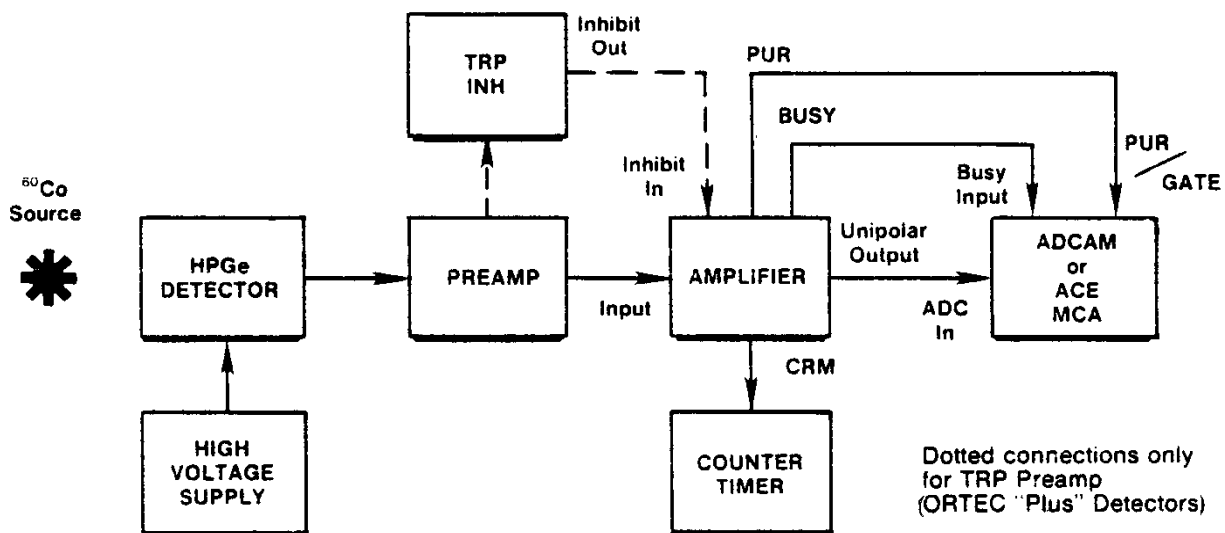


Fig. 4.11. Block Diagram for a Gamma-Ray Spectroscopy System with Pile-Up Rejection and Livetime Correction.

4.10. OPERATION WITH SEMICONDUCTOR DETECTORS

CALIBRATION OF TEST PULSER An ORTEC 419 Precision Pulse Generator, or equivalent, is easily calibrated so that the maximum pulse height dial reading (1000 divisions) is equivalent to a 10-MeV loss in a silicon radiation detector. The procedure is as follows:

- Connect the detector to be used to the spectrometer system, that is, preamplifier, main amplifier, and biased amplifier.
- Allow excitation from a source of known energy (for example, alpha particles) to fall on the detector.
- Adjust the amplifier gain and the bias level of the biased amplifier to give a suitable output pulse.
- Set the pulser Pulse Height control at the energy of the alpha particles striking the detector (e.g., set the dial at 547 divisions for a 5.47-MeV alpha particle energy).
- Turn on the pulser and use its Normalize control and attenuators to set the output due to the pulser for the same pulse height as the pulse obtained in step c. Lock the Normalize control and do not move it again until recalibration is required.

The pulser is now calibrated; the Pulse Height dial reads directly in MeV if the number of dial divisions is divided by 100.

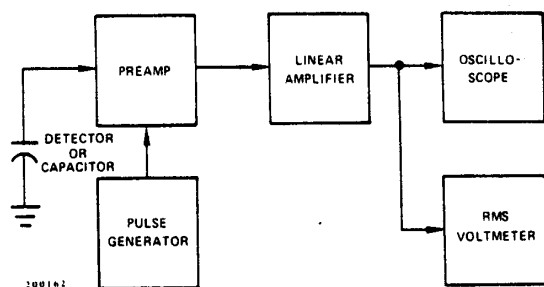


Fig. 4.12. System for Measuring Amplifier and Detector Noise Resolution.

AMPLIFIER NOISE AND RESOLUTION MEASUREMENTS

As shown in Fig. 4.12, a preamplifier, amplifier, pulse generator, oscilloscope, and wide-band rms voltmeter such as the Hewlett-Packard 3400A are required for this measurement. Connect a suitable capacitor to the input to simulate the detector capacitance desired. To obtain the resolution spread due to amplifier noise:

- Measure the rms noise voltage (E_{rms}) at the amplifier output.
- Turn on the 419 precision pulse generator and adjust the pulser output to any convenient readable voltage, E_o , as determined by the oscilloscope.

The full-width-at-half-maximum (FWHM) resolution spread due to amplifier noise is then

$$N(\text{FWHM}) = \frac{2.35 E_{rms} E_{dial}}{E_o}$$

where E_{dial} is the pulser dial reading in MeV and 2.35 is factor for rms to FWHM. For average-responding voltmeters such as the Hewlett-Packard 400D, the measured noise must be multiplied by 1.13 to calculate the rms noise.

The resolution spread will depend on the total input capacitance, since the capacitance degrades the signal-to-noise ratio much faster than the noise.

DETECTOR NOISE-RESOLUTION MEASUREMENTS

The measurement just described can be made with a biased detector instead of the external capacitor that would be used to simulate detector capacitance. The resolution spread will be larger because the detector contributes both noise and capacitance to the input. The detector noise-resolution spread can be isolated from the amplifier noise spread if the detector capacity is known, since

$$(N_{det})^2 + (N_{amp})^2 = (N_{total})^2,$$

where N_{total} is the total resolution spread and N_{amp} is the amplifier resolution spread when the detector is replaced by its equivalent capacitance.

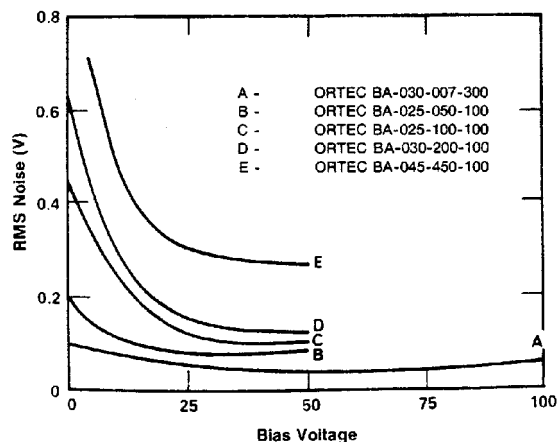


Fig. 4.13. Noise as a Function of Bias Voltage.

The detector noise tends to increase with bias voltage, but the detector capacitance decreases, thus reducing the resolution spread. The overall resolution spread will depend upon which effect is dominant. Figure 4.13 shows curves of typical noise-resolution spread versus bias voltage, using data from several ORTEC silicon surface-barrier semiconductor radiation detectors.

AMPLIFIER NOISE-RESOLUTION MEASUREMENTS USING MCA Probably the most convenient method of making resolution measurements is with a pulse height analyzer as shown by the setup illustrated in Fig. 4.14.

The amplifier noise-resolution spread can be measured directly with a pulse height analyzer and the mercury pulser as follows:

- Select the energy of interest with an ORTEC 419 Precision Pulse Generator. Set the amplifier and biased amplifier gain and bias level controls so that the energy is in a convenient channel of the analyzer.
- Calibrate the analyzer in keV per channel, using the pulser; full scale on the pulser dial is 10 MeV when calibrated as described above.
- Obtain the amplifier noise-resolution spread by measuring the FWHM of the pulser peak in the spectrum.

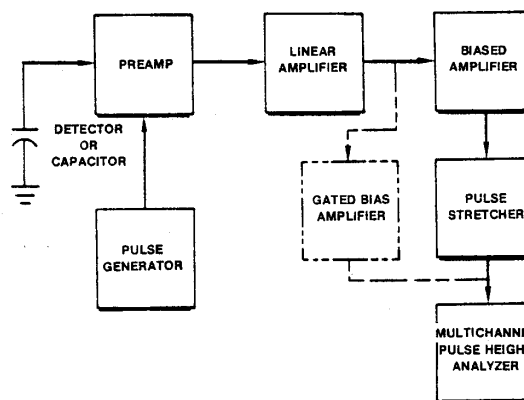


Fig. 4.14. System for Measuring Resolution with a Pulse Height Analyzer.

The detector noise-resolution spread for a given detector bias can be determined in the same manner by connecting a detector to the preamplifier input. The amplifier noise resolution spread must be subtracted as described in Section 4.10, "Detector Noise-Resolution Measurements." The detector noise will vary with detector size and bias conditions and possibly with ambient conditions.

CURRENT-VOLTAGE MEASUREMENTS FOR Si AND Ge DETECTORS The amplifier system is not directly involved in semiconductor detector current-voltage measurements, but the amplifier serves to permit noise monitoring during the setup. The detector noise measurement is a more sensitive method of determining the maximum detector voltage than a current measurement and should be used because the noise increases more rapidly than the reverse current at the onset of detector breakdown. Make this measurement in the absence of a source.

Figure 4.15 shows the setup required for current-voltage measurements. An ORTEC 428 Bias Supply is used as the voltage source. Bias voltage should be applied slowly and reduced when noise increases rapidly as a function of applied bias. Figure 4.16 shows several typical current-voltage curves for ORTEC silicon surface-barrier detectors.

When it is possible to float the microammeter at the detector bias voltage, the method of detector current measurement shown by the dashed lines in

Fig. 4.15 is preferable. The detector is grounded as in normal operation, and the microammeter is connected to the current monitoring jack on the 428 detector bias supply.

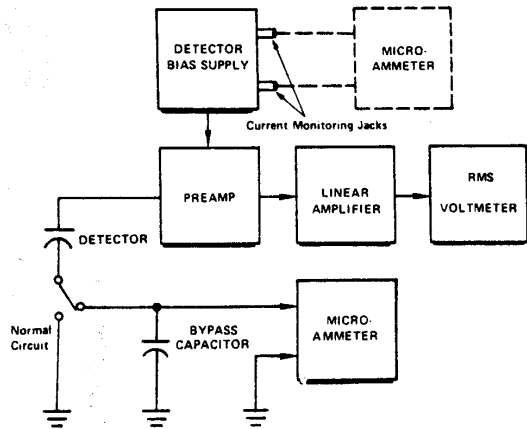


Fig. 4.15. System for Detector Current and Voltage Measurements.

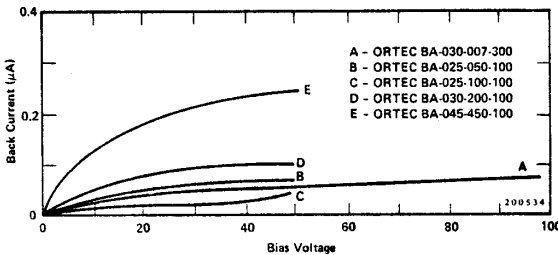


Fig. 4.16. Silicon Detector Back Current vs Bias Voltage.

4.11. OPERATION IN SPECTROSCOPY SYSTEMS

HIGH-RESOLUTION ALPHA-PARTICLE SPECTROSCOPY SYSTEM The block diagram of a high-resolution spectroscopy system for measuring natural alpha particle radiation is shown in Fig. 4.17. Since natural alpha radiation occurs only above several MeV, an ORTEC 444 Biased Amplifier is used to suppress the unused portion of the spectrum; the same result can be obtained by using digital suppression on the MCA in many cases. Alpha-particle resolution is obtained in the following manner:

- a. Use appropriate amplifier gain and minimum biased amplifier gain and bias level. Accumulate the alpha peak in the MCA.
- b. Slowly increase the bias level and biased amplifier gain until the alpha peak is spread over 5 to 10 channels and the minimum- to maximum-energy range desired corresponds to the first and last channels of the MCA.
- c. Calibrate the analyzer in keV per channel using the pulser and the known energy of the alpha peak (see Section 4.10, "Calibration of Test Pulser") or two known energy alpha peaks.
- d. Calculate the resolution by measuring the number of channels at the FWHM level in the peak and converting this to keV.

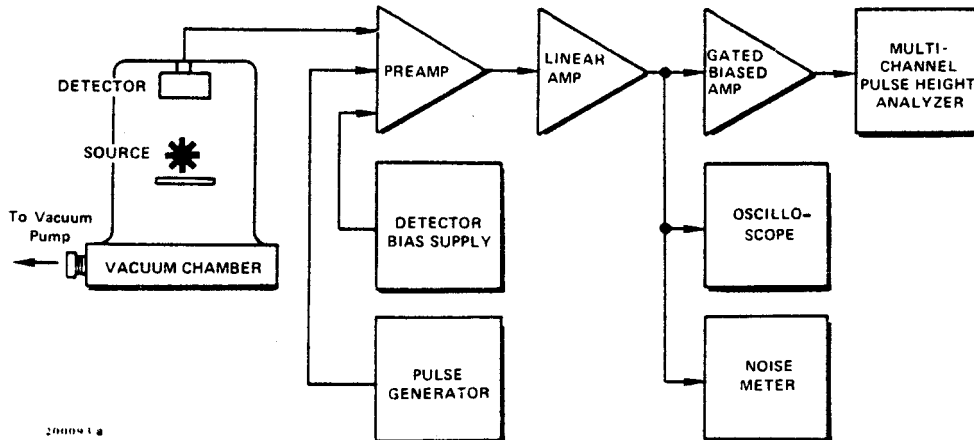


Fig. 4.17. System for High-Resolution Alpha-Particle Spectroscopy.

HIGH-RESOLUTION GAMMA-RAY SPECTROSCOPY SYSTEM

A high-resolution gamma-ray spectroscopy system block diagram is shown in Fig. 4.18. Although a biased amplifier is not shown (an analyzer with more channels being preferred), it can be used if the only analyzer available has fewer channels and only higher energies are of interest.

When germanium detectors that are cooled by a liquid nitrogen cryostat are used, it is possible to obtain resolutions from about 1 keV FWHM up to 4 keV (depending on the energy of the incident radiation and the size and quality of the detector). Reasonable care is required to obtain such results. Some guidelines for obtaining optimum resolution are:

- Keep interconnection capacities between the detector and preamplifier to an absolute minimum (no long cables).
- Keep humidity low near the detector-preamplifier junction.
- Operate the amplifier with the shaping time that provides the best signal-to-noise ratio.
- Operate at the highest allowable detector bias to keep the input capacity low.

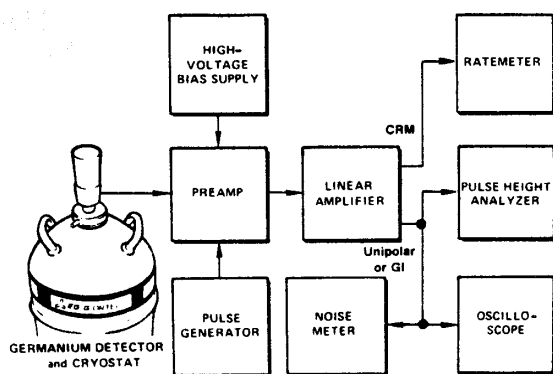


Fig. 4.18. System for High-Resolution Gamma Spectroscopy.

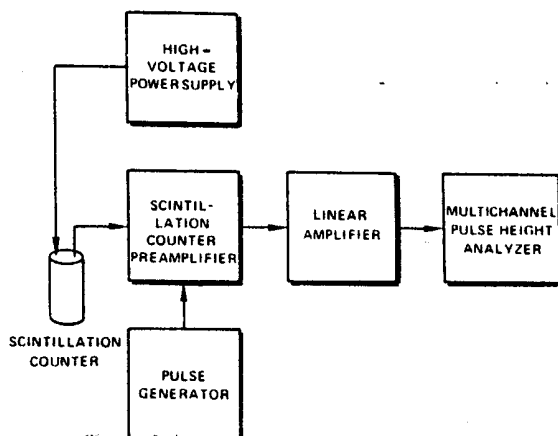


Fig. 4.19. Scintillation-Counter Gamma Spectroscopy System.

SCINTILLATION-COUNTER GAMMA SPECTROSCOPY SYSTEMS

The ORTEC 672 can be used in scintillation-counter spectroscopy systems as shown in Fig. 4.19. The amplifier shaping time constants should be selected in the region of 0.5 to 1 μ s for NaI or plastic scintillators. For scintillators having longer decay times, longer time constants should be selected.

X-RAY SPECTROSCOPY USING PROPORTIONAL COUNTERS

Space charge effects in proportional counters, operated at high gas amplification, tend to degrade the resolution capabilities drastically at x-ray energies, even at relatively low counting rates. By using a high-gain low-noise amplifying system and lower gas amplification, these effects can be reduced and a considerable improvement in resolution can be obtained. The block diagram in Fig. 4.20 shows a system of this type. Analysis can be accomplished by simultaneous acquisition of all data on a multichannel analyzer or counting a region of interest in a single-channel analyzer window with a counter and timer or counting ratemeter.

4.12. OTHER EXPERIMENTS

Block diagrams illustrating how the 672 and other ORTEC modules can be used for experimental setups for various other applications are shown in Figs. 4.21 through 4.24.

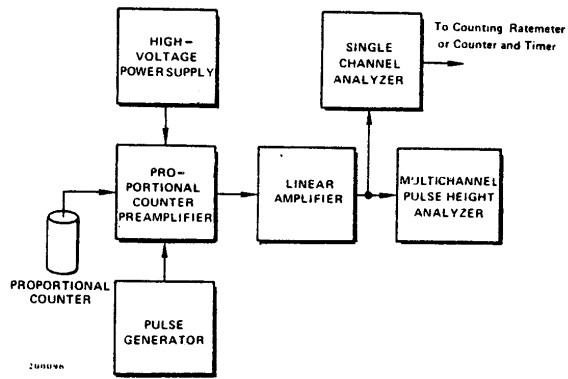


Fig. 4.20. High-Resolution X-Ray Energy Analysis System Using a Proportional Counter.

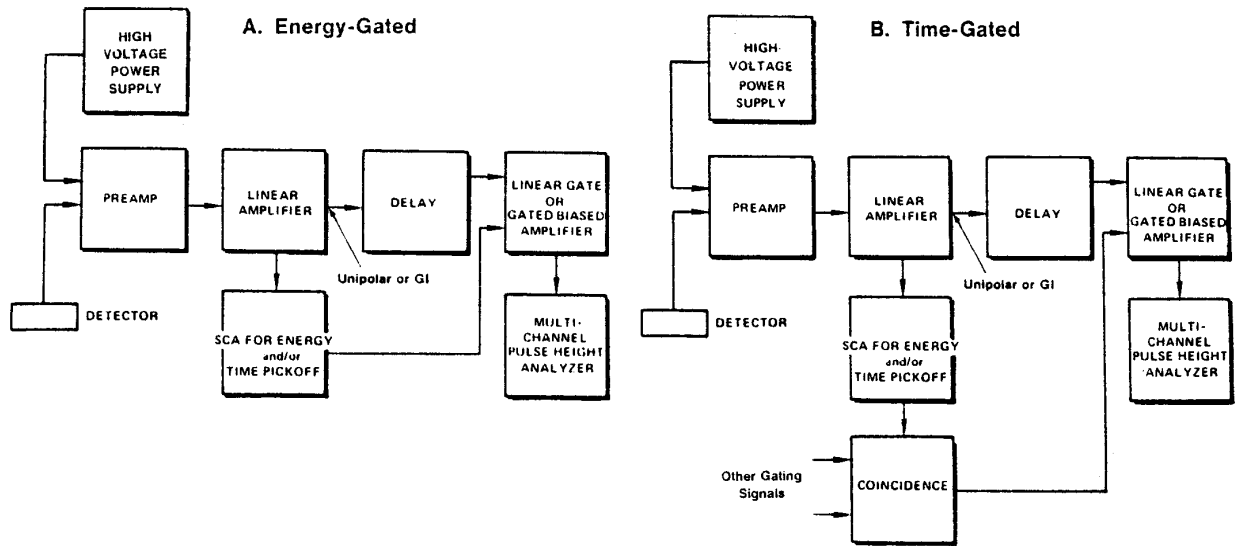


Fig. 4.21. General System Arrangement for Gating Control.

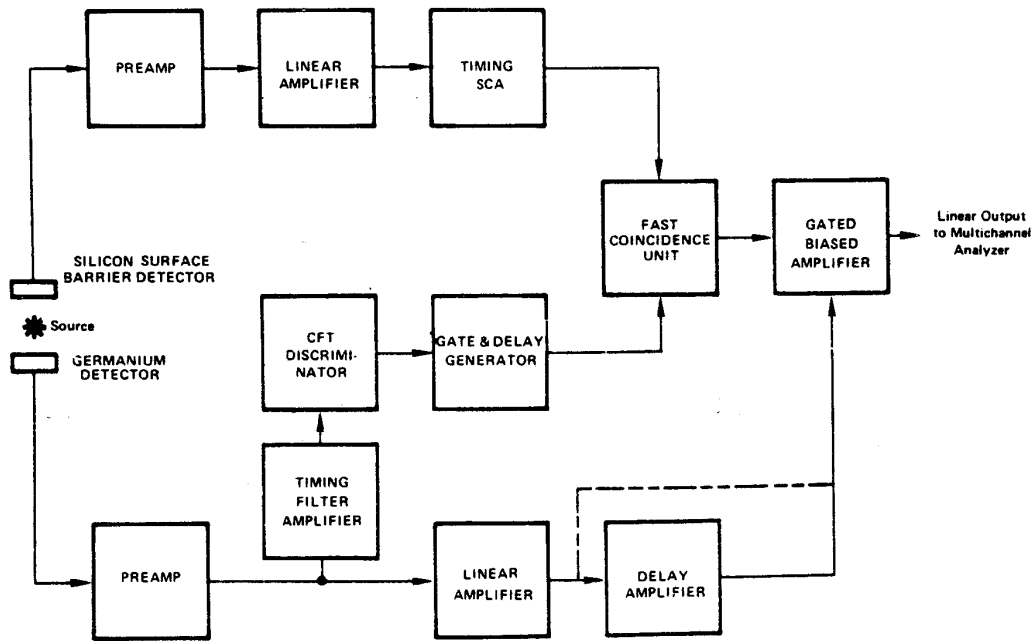


Fig. 4.22. Gamma-Ray Charged-Particle Coincidence Experiment.

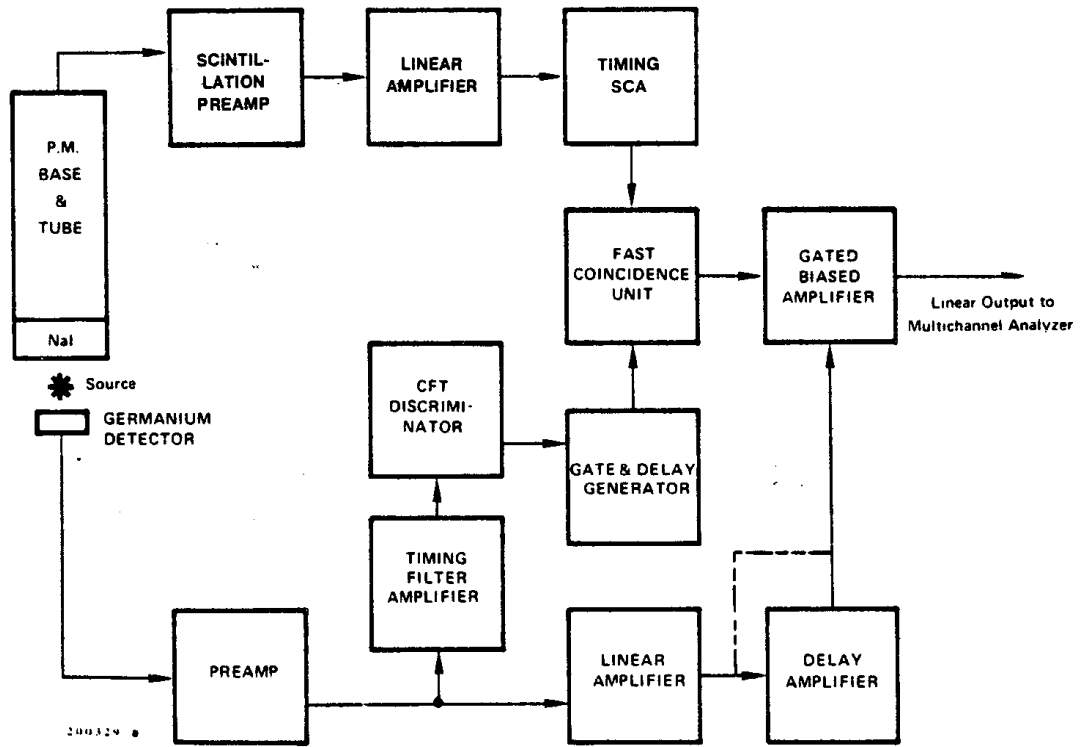


Fig. 4.24. Gamma-Gamma Coincidence Experiment.

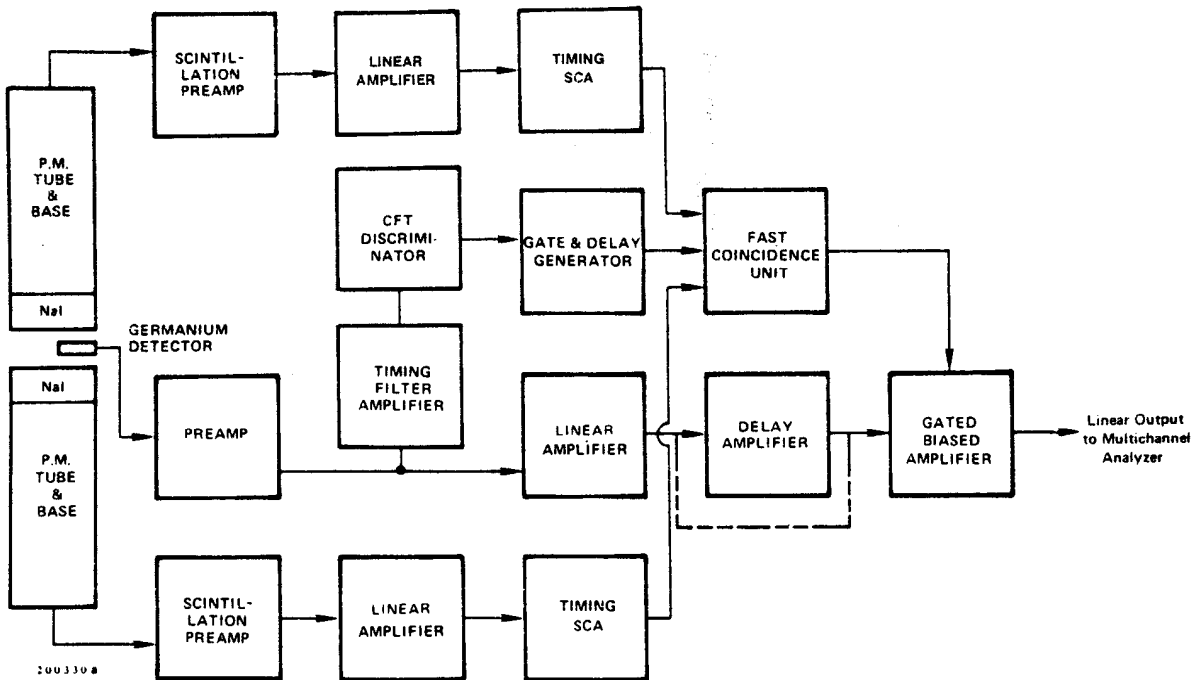


Fig. 4.23. Gamma-Ray Pair Spectroscopy.

5. MAINTENANCE

5.1. TEST EQUIPMENT REQUIRED

The following test equipment should be utilized to adequately test the specifications of the 672 Spectroscopy Amplifier:

1. ORTEC 419 Precision Pulse Generator or 448 Research Pulser.
2. Tektronix 465, 475, or 485 Series Oscilloscope or equivalent.
3. Hewlett-Packard 3400A RMS Voltmeter.

5.2. PULSER TEST⁶

Coarse Gain	1K
Gain	1.5
Input Polarity	Positive
Shaping Time Constant	2 μ s
BLR	PZ
UNI Shaping	Gaussian
PZ	MAN

- a. Connect a positive pulser output to the 672 input and adjust the pulser to obtain +10 V at the 672 Unipolar output. This should require an input pulse of 6.6 mV, using a 100 Ω terminator at the input. Switch Unipolar Mode to Triangle. This should also be 10 V.
- b. Measure the positive lobe of the Bipolar output. This should also be +10 V.
- c. Change the Input polarity switch to Neg and then back to Pos while monitoring the outputs for a polarity inversion. The negative output should clamp at -1 V.
- d. Auto PZ Check - Adjust a 448 pulser or equivalent pulser with the decay time set for 50 μ s

to obtain a positive 9-V pulse at the Unipolar output. A 419 pulser with its longer decay time can also be used, but the time needed for the Auto PZ process may be quite long at the 60- or 70-Hz repetition rate. Select Auto PZ and push the button. The red LED will go on and then turn off when the unipolar output pulse has been pole-zero compensated. The quality of the PZ can be observed by increasing the scope input sensitivity to 10 or 20 mV/cm with the LIMIT button pushed. The LIMIT push-button switch inserts a diode limiter in series with the front-panel Unipolar Output connector to prevent overloading the input of the oscilloscope.

Change the Shaping Time switch to 0.5 through 10 μ s. At each setting the PZ should be close to being correct. A push of the Auto PZ button should complete the pole-zero compensation.

- e. Decrease the Coarse Gain switch stepwise from 1K to 5 and ensure that the output amplitude changes by the appropriate amount for each step. Return the Coarse Gain switch to 1 K.
- f. Decrease the Fine Gain control from 1.5 to 0.5 and check to see that the output amplitude decreases by a factor of 3. Return the Fine Gain control to maximum at 1.5.
- g. With the Shaping Time switch set for 1 μ s, measure the time to the peak on the unipolar output pulse; this should be 2.2T.
- h. Change the Shaping Time switch to 0.25 through 6 μ s. At each setting, check to see that the time to the unipolar peak is 2.2T. Return the switch to 1 μ s.

OVERLOAD TESTS Start with maximum gain, $\tau=2\mu$ s, and a +10 V output amplitude. Increase the pulser output amplitude by X1000 and observe that the unipolar output returns to within 200 mV of the baseline within 27 μ s after the application of a single pulse from the pulser. It will probably be necessary to vary the PZ Adj control on the front panel in order to cancel the pulser pole and minimize the time required for return to the baseline.

⁶See IEEE Standards, No. 301-1976.

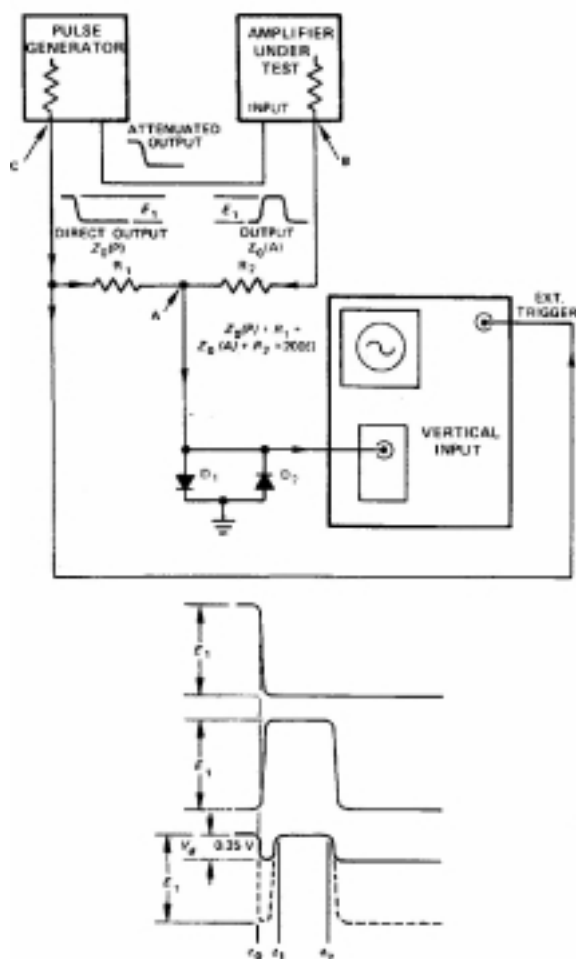


Fig. 5.1. Circuit Used to Measure Nonlinearity.

LINEARITY The integral nonlinearity of the 672 can be measured by the technique shown in Fig. 6.1. In effect, the negative pulser output is subtracted from the positive amplifier output to cause a null point that can be measured with excellent sensitivity. The pulser output must be varied between 0 and 10V, which usually requires an external control source for the pulser. The amplifier gain and the pulser attenuator must be adjusted to measure 0 V at the null point when the pulser output is 10 V. The variation in the null point as the pulser is reduced gradually from 10 V to 0 V is a measure of the nonlinearity. Since the subtraction network also acts as a voltage divider, this variation must be less than

$$(10 \text{ V full scale}) \times (\pm 0.025\% \text{ maximum nonlinearity}) \\ \times (1/2 \text{ for divider network}) = \pm 1.25 \text{ V}$$

for the maximum null-point variation.

OUTPUT LOADING Use the test setup of Fig. 5.1. Adjust the amplifier output to 10 V and observe the null point when the front panel output is terminated in 100Ω. The change should be <5 mV.

NOISE Measure the noise at the amplifier Unipolar output with maximum amplifier gain and 2-μs shaping time. Using a true rms voltmeter, the noise should be less than 5 μV x 1500 (gain), or 7.5 mV.

For an average responding voltmeter, the noise reading would have to be multiplied by 1.13 to calculate the rms noise. The input must be terminated in 100Ω during the noise measurements

5.3. SUGGESTIONS FOR TROUBLESHOOTING

In situations where the 672 is suspected of a malfunction, it is essential to verify such malfunction in terms of simple pulse generator impulses at the input. The 672 must be disconnected from its position in any system, and routine diagnostic analysis performed with a test pulse generator and an oscilloscope. It is imperative that testing not be performed with a source and detector until the amplifier performs satisfactorily with the test pulse generator.

The testing instructions in Section 5.2 should provide assistance in locating the region of trouble and repairing the malfunction. The two side plates can be completely removed from the module to enable oscilloscope and voltmeter observations.

5.4. FACTORY REPAIR

This instrument can be returned to the ORTEC factory for service and repair at a nominal cost. Our standard procedure for repair ensures the same quality control and checkout that are used for a new instrument. Always contact Customer Services at ORTEC, (865) 483-2231, before sending in an instrument for repair to obtain shipping instructions and so that the required Return Authorization Number can be assigned to the unit. This number should be marked on the address label and on the package to ensure prompt attention when the unit reaches the factory.

5.5. TABULATED TEST POINT VOLTAGES

The voltages given in Table 5.1 are intended to indicate typical dc levels that can be measured on the printed circuit board. In some cases the circuit will perform satisfactorily even though, due to component tolerances, there may be some voltage measurements that differ slightly from the listed values. Therefore the tabulated values should not be interpreted as absolute voltages but are intended to serve as an aid in troubleshooting.

Table 5.1. Typical dc Voltages*

Location	Voltage
TP1	±25 mV
TP2	±25 mV
TP3	±25 mV
TP4	±10 mV
TP5	±10 mV
TP6	±25 mV
TP7	±40 mV
TP8	±50 mV
TP9	-12.5 V
TP10	-11.4 V
TP11	+0.05 to +0.2 V
TP12	-0.5 V
TP13	HC Logic 0
TP14	HC Logic 1
TP15	±25 mV
TP16	-0.4 V
TP17	+0.05 to +0.2 V
TP18	-11.3 V
UNI Out	±5 mV
BI Out	±10 mV

*All voltages measured with no input signal, with the input terminated in 100 Ω, and all controls set fully clockwise at maximum.

**Bin/Module Connector Pin
Assignments For Standard
Nuclear Instrument Modules per
DOE/ER-0457T.**

Pin	Function	Pin	Function
1	+3 V	23	Reserved
2	- 3 V	24	Reserved
3	Spare bus	25	Reserved
4	Reserved bus	26	Spare
5	Coaxial	27	Spare
6	Coaxial	*28	+24 V
7	Coaxial	*29	- 24 V
8	200 V dc	30	Spare bus
9	Spare	31	Spare
*10	+6 V	32	Spare
*11	- 6 V	*33	117 V ac (hot)
12	Reserved bus	*34	Power return ground
13	Spare	35	Reset (Scaler)
14	Spare	36	Gate
15	Reserved	37	Reset (Auxiliary)
*16	+12 V	38	Coaxial
*17	- 12 V	39	Coaxial
18	Spare bus	40	Coaxial
19	Reserved bus	*41	117 V ac (neutral)
20	Spare	*42	High-quality ground
21	Spare	G	Ground guide pin
22	Reserved		

Pins marked (*) are installed and wired in ORTEC's 4001A and 4001C Modular System Bins.